Edwards Fogarty[®] Catheters

The only name in clot management backed by over 50 years of worldwide experience.



For over 50 years, clinicians have cou to address the real problem of clot

The reason is simple. Through collaboration and clinical learning, Edwards Lifesciences has continued to evolve and advance its Fogarty[°] catheters to meet the exacting needs of specialists for vascular procedures. Their proven performance and manufacturing quality has led surgeons to specify genuine Fogarty[°] catheters – available exclusively from Edwards – for millions of procedures worldwide.

Fresh, soft emboli/thrombi

Fogarty[®] arterial embolectomy catheter

Removal of fresh, soft emboli and thrombi from vessels in the arteries of the non-central circulatory system

- Hand-tied using a recessed winding technique to secure balloon to catheter
- Balloon exhibits a symmetry that exerts uniform contact with vessel walls for even pressure and precise traction
- Rounded tip promotes easy insertion for reduced trauma

Fogarty° Fortis arterial embolectomy catheter

Removal of fresh, soft emboli and thrombi from vessels in the arteries of the non-central circulatory system

- Can withstand greater pull force than the Fogarty^{*} arterial embolectomy catheter
- The addition of barium to the catheter tubing provides radiopacity

Fresh, soft emboli/thrombi, temporary occlusion, infusion

Fogarty[®] thru-lumen embolectomy catheter

Multi-use guidewire-compatible design enables removal of fresh, soft emboli and thrombi from vessels in the arterial system as well as temporary occlusion of blood vessels, infusion of fluids, and blood sampling

- Fogarty[®] balloon conforms to the vessel to maximize clot removal while minimizing trauma
- Provides the solution for distal thrombus removal in the arterial system
- The addition of barium to the catheter shaft material and stainless steel markers proximal and distal to the balloon provide radiopacity*

* 3F Thru-Lumen embolectomy catheter body does not include barium







nted on real Fogarty[®] catheters management.

Mature, adherent thrombi

Fogarty[®] Corkscrew catheter

Indicated for removal of emboli and thrombi from either native vessels or synthetic grafts in the arterial system

- For extracting thrombus that is too resistant for an elastomeric balloon
- Corkscrew design increases the surface area for entrapping fibrous material
- Radiopaque and flexible

Fogarty[®] graft thrombectomy catheter

Designed for removing organized thrombus and adherent thrombotic material from synthetic bypass grafts, specifically, AV fistulas for dialysis access and grafts for iliofemoral bypass


- For thrombus that is too resistant for an elastomeric balloon
- Provides a larger diameter range than the Fogarty[®] Corkscrew catheter
- The pulling force can be varied by adjusting the diameter

Venous thrombi

Fogarty[®] venous thrombectomy catheter

Designed for removing thrombi in the venous system

• The long, soft tip permits advancement past the venous valves without undue trauma

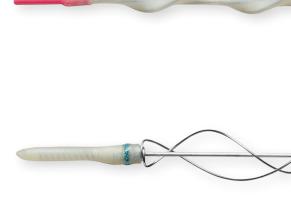



You'll want a real



Fogarty® clot management catheters										
Model	Length (cm)	Catheter size (F)	Membrane/spiral diameter max. extended (mm)	Membrane/spiral diameter max. retracted (mm)						
Fogarty [®] Corkscrew catheter										
140806	80	4	2.7	6						
140808	80	5	3.1	8						
1408010	80	6	3.3	10						
Fogarty' graft thrombectomy catheters										
160245F	50	5	5.0	16						
160246F	50	6	6.0	18						

Model	Length (cm)	Catheter size (F)	Inflated balloon diameter (mm)	Max. size of deflated balloon (F)	Max. liquid capacity (ml)	Max. recom- mended pull force (lbs.)	Compatible guidewire (in.)			
Fogarty [®] arterial embolectomy catheters (pouches, 5 pack)										
120602FP	60	2	4	3.4	-	0.5	-			
120403FP	40	3	5	4.3	0.20	0.7	-			
120803FP	80	3	5	4.3	0.20	0.7	-			
120404FP	40	4	9	5.0	0.75	1.5	-			
120804FP	80	4	9	5.0	0.75	1.5	-			
120805FP	80	5	11	6.0	1.50	2.0	-			
120806FP	80	6	13	7.0	2.00	2.5	-			
120807FP	80	7	14	8.0	2.50	3.5	-			
Fogarty [®] Fortis arterial embolectomy catheters										
120404FF	40	4	9	5.0	0.75	-	-			
120804FF	80	4	9	5.0	0.75	-	-			
Fogarty [®] thru-lumen embolectomy catheters										
12TLW403F	40	3	5	4.5	0.15	0.7	0.018			
12TLW803F	80	3	5	4.5	0.15	0.7	0.018			
12TLW1203F	120	3	5	4.5	0.15	0.7	0.018			
12TLW404F	40	4	9	5.0	0.50	1.5	0.025			
12TLW804F	80	4	9	5.0	0.50	1.5	0.025			
12TLW405F35	40	5.5	11	6.0	0.90	2.0	0.035			
12TLW805F35	80	5.5	11	6.0	0.90	2.0	0.035			
12TLW406F	40	6	13	7.0	1.40	2.5	0.035			
12TLW806F	80	6	13	7.0	1.40	2.5	0.035			
12TLW807F	80	7	14	8.0	1.60	3.5	0.038			
Fogarty [®] venous thrombectomy catheters										
320806F	80	6	12	7.7	1.50	-	-			
320808F	80	8	13	8.8	2.25	-	-			
32080810F	80	8	19	10.0	4.00	_	-			

Genuine Fogarty[®] catheters are available exclusively from Edwards Lifesciences. To learn more, visit edwards.com/fogarty or call 1-800-FOGARTY.

CAUTION: Federal (United States) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, precautions and adverse events.

Edwards, Edwards Lifesciences, the stylized E logo, Fogarty, Fogarty Corkscrew, Fogarty Fortis, and 1-800-FOGARTY are trademarks or service marks of Edwards Lifesciences Corporation or its affiliates. All other trademarks are the property of their respective owners.

© 2020 Edwards Lifesciences Corporation. All rights reserved. PP--US-0590 v3.0

Edwards Lifesciences • One Edwards Way, Irvine CA 92614 USA • edwards.com

