


# Cordis **STORQ**<sup>TM</sup>

Steerable Guidewires



Steerability and Torque,  
experience the Control

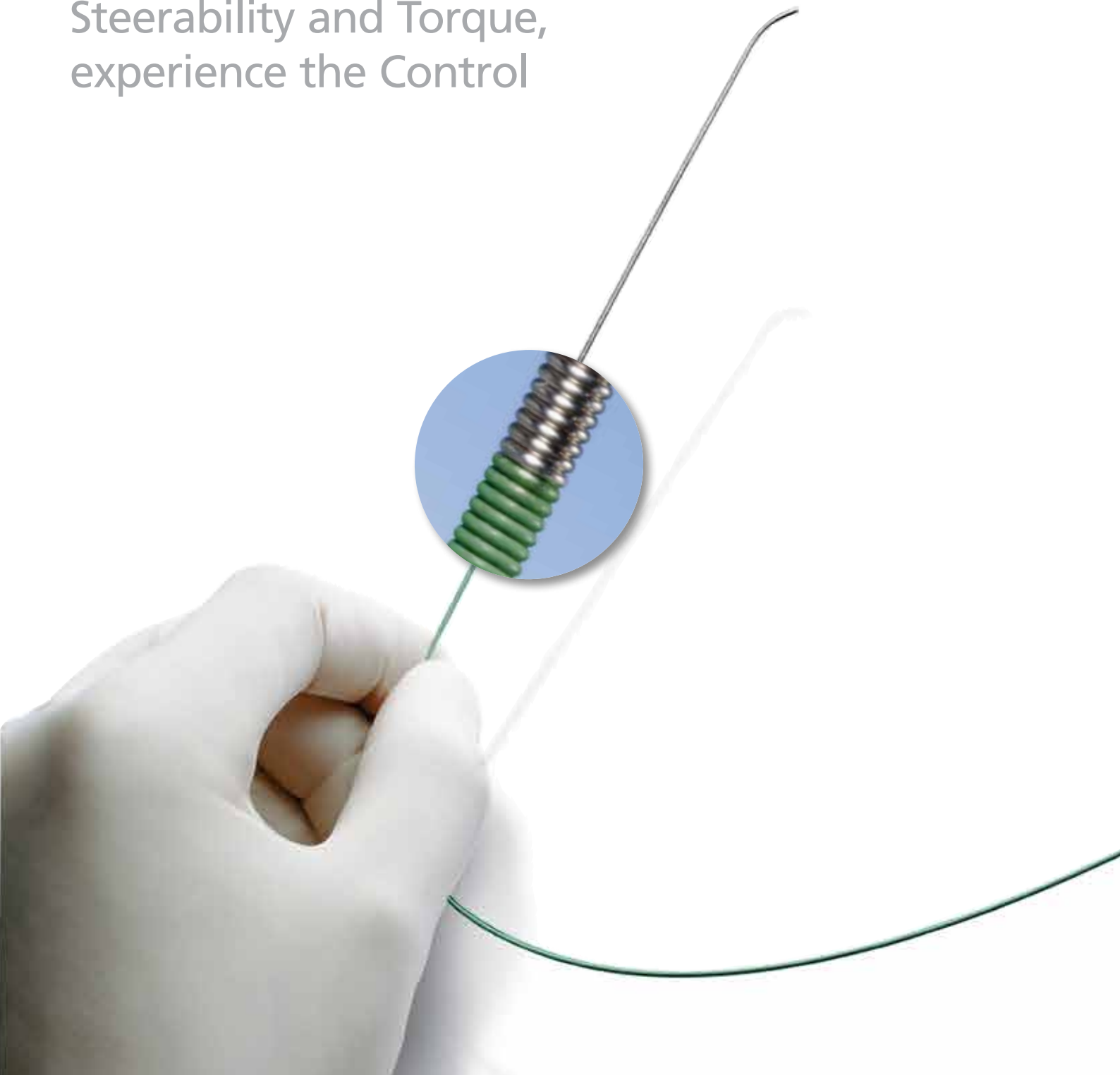

# Cordis **STORQ**<sup>TM</sup>

Steerable Guidewires

## The **STORQ**<sup>TM</sup> Guidewire brings you a softer tip and sturdy support

### Steerability and Torque, Experience the control

- Designed STORQ<sup>TM</sup> guidewire's atraumatic tip segment is soft and flexible to guard against intimal damage or plaque disruption during device delivery and catheter exchange.
- Length of core taper provides three different distal flexibilities - soft, supersoft, and standard - for accurate placement in diverse anatomy.
- Preshaped angled, modified J, and straight tip shapes provide physician with procedural flexibility.



### Excellent Tip and Body Visualisation


- Radiopaque distal tip and radiolucent shaft permit clear fluoroscopic visualisation that allows precise placement.




## Different Taper Lengths for Procedural Flexibility

If you want...	Then select...	Taper length
<b>Firm support</b>	Standard	10 cm
<b>Multifunctional capability</b>	Soft	13 cm
<b>Maximum flexibility</b>	SuperSoft	25 cm

PTFE and SLX™  
Guidewire Coatings

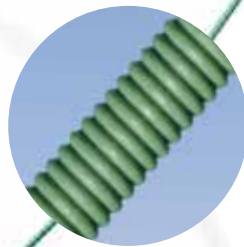

Torque Device

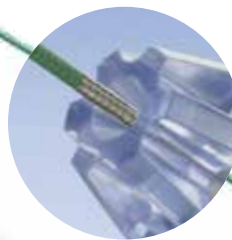

## PTFE and SLX™ Coatings for Superb Lubricity

- PTFE (polytetrafluoroethylene) and SLX™ (silicone) coatings cover most of the guidewire, creating a lubricious surface for improved catheter tracking and guidewire movement.
  - Proximal segment is uncoated for secure attachment of torque device.
  - PTFE coating ends 10 cm from distal tip, while SLX™ coating continues to tip to facilitate crossing tight or tortuous lesions.
- Designed with durable PTFE and SLX™ coating.

## Excellent Support

- STORQ™ provides enhanced shaft support and stability for stent and balloon catheter delivery.
- STORQ™ maintains excellent torqueability and steerability for controlled access to target lesions.
- Design features:
  - Rigid, one-piece, stainless steel core for excellent torque response to prevent whipping.
  - Single outer coil that runs entire length of wire for more flexibility.

# Cordis **STORQ**<sup>TM</sup>

Steerable Guidewires

Packaged five per box  
Torque device included with each guidewire



## Ordering information

Diameter (inch)	Tip Flexibility	Tip Shape	Length (cm)	Catalogue Number
.035	SuperSoft	Straight	180	503-256
.035	SuperSoft	Angled	180	503-256J
.035	Soft	Straight	180	503-356
.035	Soft	Angled	180	503-356J
.035	Soft	Modified J	180	503-356MJ
.035	Soft	Straight	300	503-356X
.035	Soft	Angled	300	503-356Y
.035	Soft	Modified J	300	503-356MY
.035	Standard	Straight	180	503-456
.035	Standard	Angled	180	503-456J
.035	Standard	Modified J	180	503-456MJ
.035	Standard	Straight	300	503-456X
.035	Standard	Angled	300	503-456Y
.035	Standard	Modified J	300	503-456MY