


Single-Loop & Triple-Loop Retrieval Kits




Redefining your retrieval options

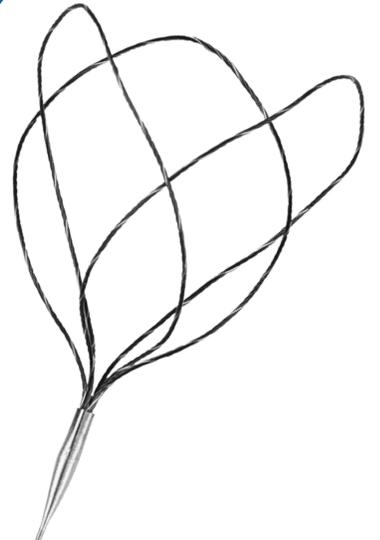

Nitinol wire
for 1:1 torque
control and kink
resistance

**Three
independent
loops**
for complete
vessel coverage

**Platinum strands
or gold-plated
tungsten**
for visibility
when snaring
retrieval hook

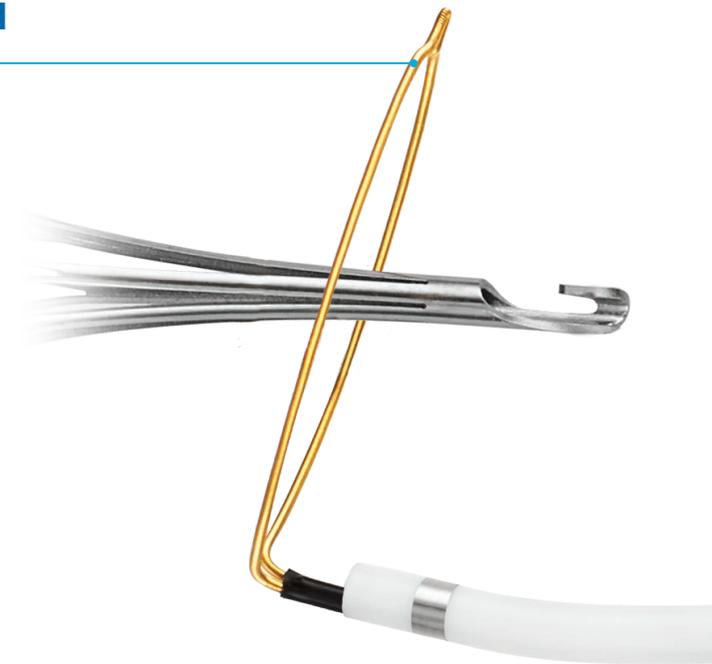

**Holes around
the dilator**
for even distribution
of contrast during
power injection

**15° curve at
the distal tip**
for directional
control




The only two kits designed for removal of all IVC filters retrievable from a jugular approach¹

With a rigid sheath, enhanced visualization, hemostatic control and two snare options, Argon's Single-Loop and Triple-Loop Retrieval Kits are designed with the rigor of the IVC filter retrieval procedure in mind.

Radiopaque bands

increase visibility for control when capturing and sheathing over the IVC filter

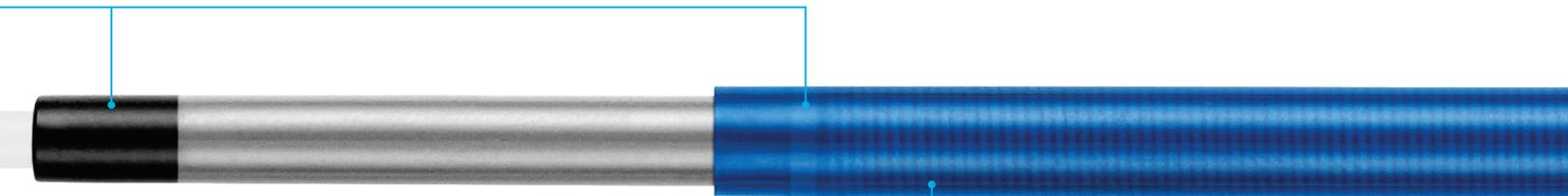

Triple-Loop Retrieval Kit

Coil-reinforced outer sheath

for strength to detach filter legs from the IVC wall



Single-Loop Retrieval Kit


Ordering Information

Single-Loop Retrieval Kit			
CATALOG #	INNER SHEATH	OUTER SHEATH	SNARE CATHETER
380911120	9F (ID) x 64cm	11F (ID) x 60cm	7F (OD) x 76cm

Triple-Loop Retrieval Kit			
CATALOG #	INNER SHEATH	OUTER SHEATH	SNARE CATHETER
380911330	9F (ID) x 64cm	11F (ID) x 60cm	7F (OD) x 76cm

Each Kit Includes

- A. 11F (ID) x 60cm Outer Sheath
- B. 9F (ID) x 64cm Inner Sheath
- C. 8F (OD) x 75cm Dilator
- D. Hemostasis Valve with Sideport
- E. 3-Way High Pressure Stopcock
- F. 7F (OD) x 76cm Snare Catheter with Tuohy-Borst Y-Port Adapter, 30mm x 93cm Triple-Loop Snare with Torque Handle
- OR
- G. 7F (OD) x 76cm Snare Catheter with Tuohy-Borst Y-Port Adapter, 20mm x 93cm Single-Loop Snare with Torque Handle

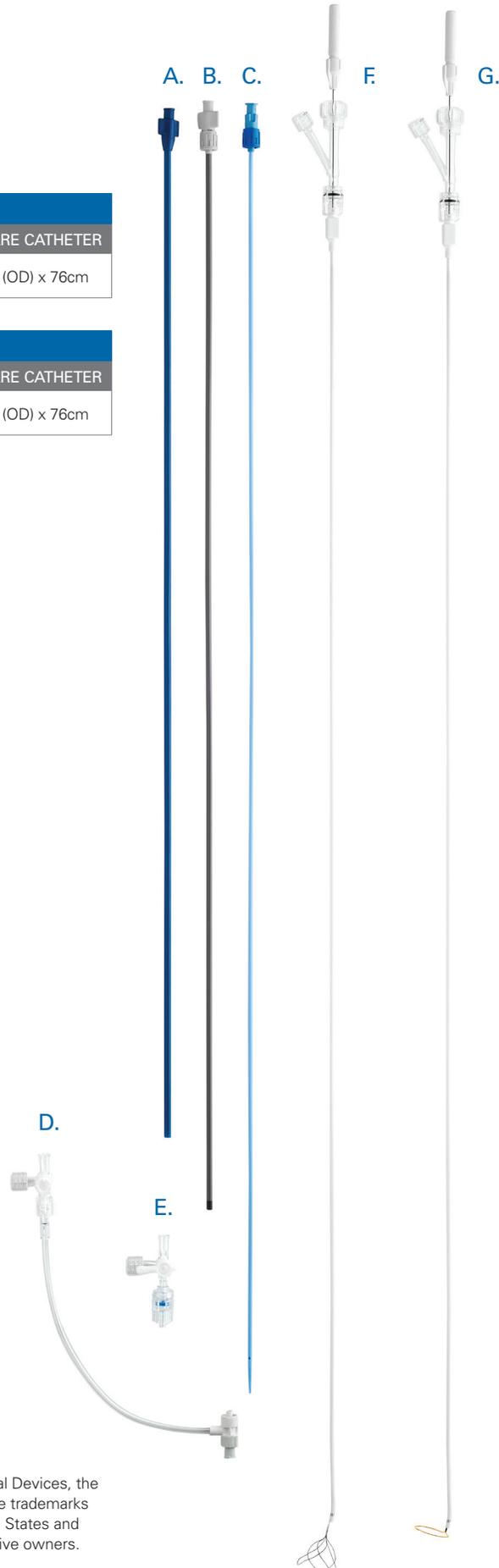

1. K070514, K152560, K073374, K112185

© 2019 Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries. All other trademarks are the property of their respective owners.

2600 Dallas Parkway, Suite 440, Frisco, TX 75034
 P: 800.927.4669 F: 903.677.9396 argonmedical.com

ARGON
 MEDICAL DEVICES

96-6500-09 1219A