



Single-Loop & Triple-Loop Retrieval Kits

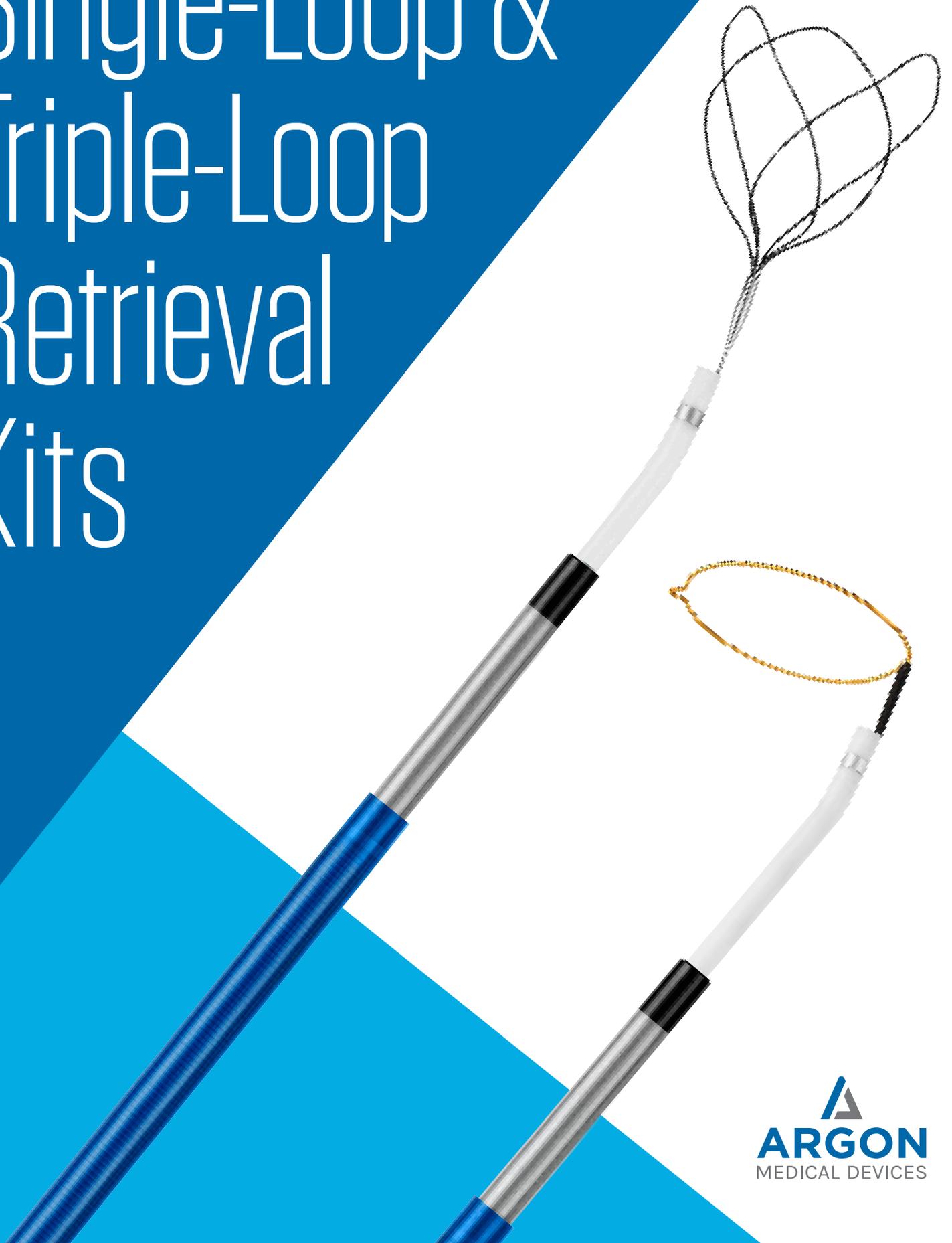

Redefining your retrieval options

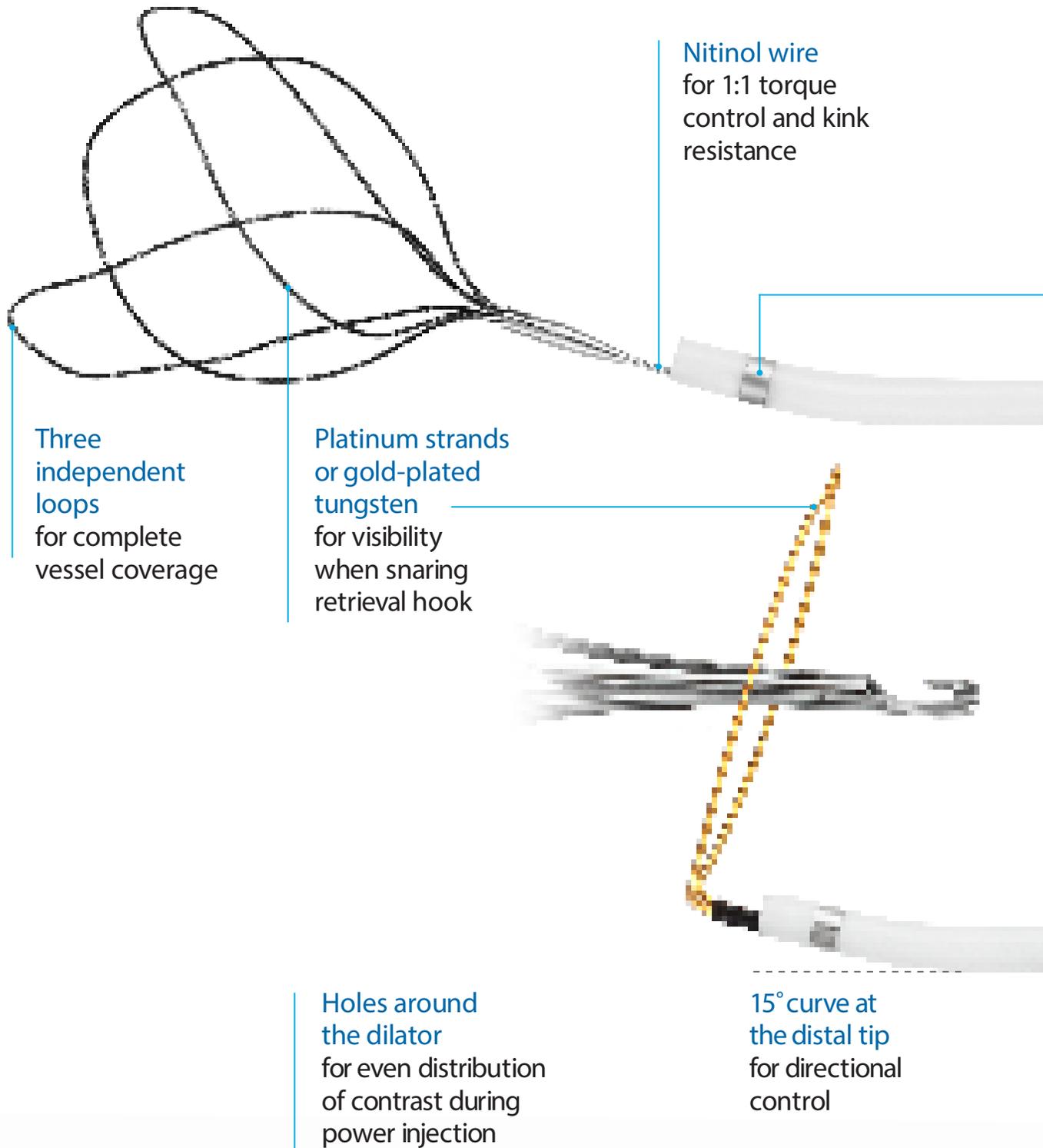


The only single-loop and triple-loop kits designed for removal of all IVC filters retrievable from a jugular approach.¹

With a rigid sheath, enhanced visualization, hemostatic control and two snare options, Argon's Single-Loop and Triple-Loop Retrieval Kits are designed with the rigor of

Radiopaque bands

increase visibility for control when capturing and sheathing over the IVC filter

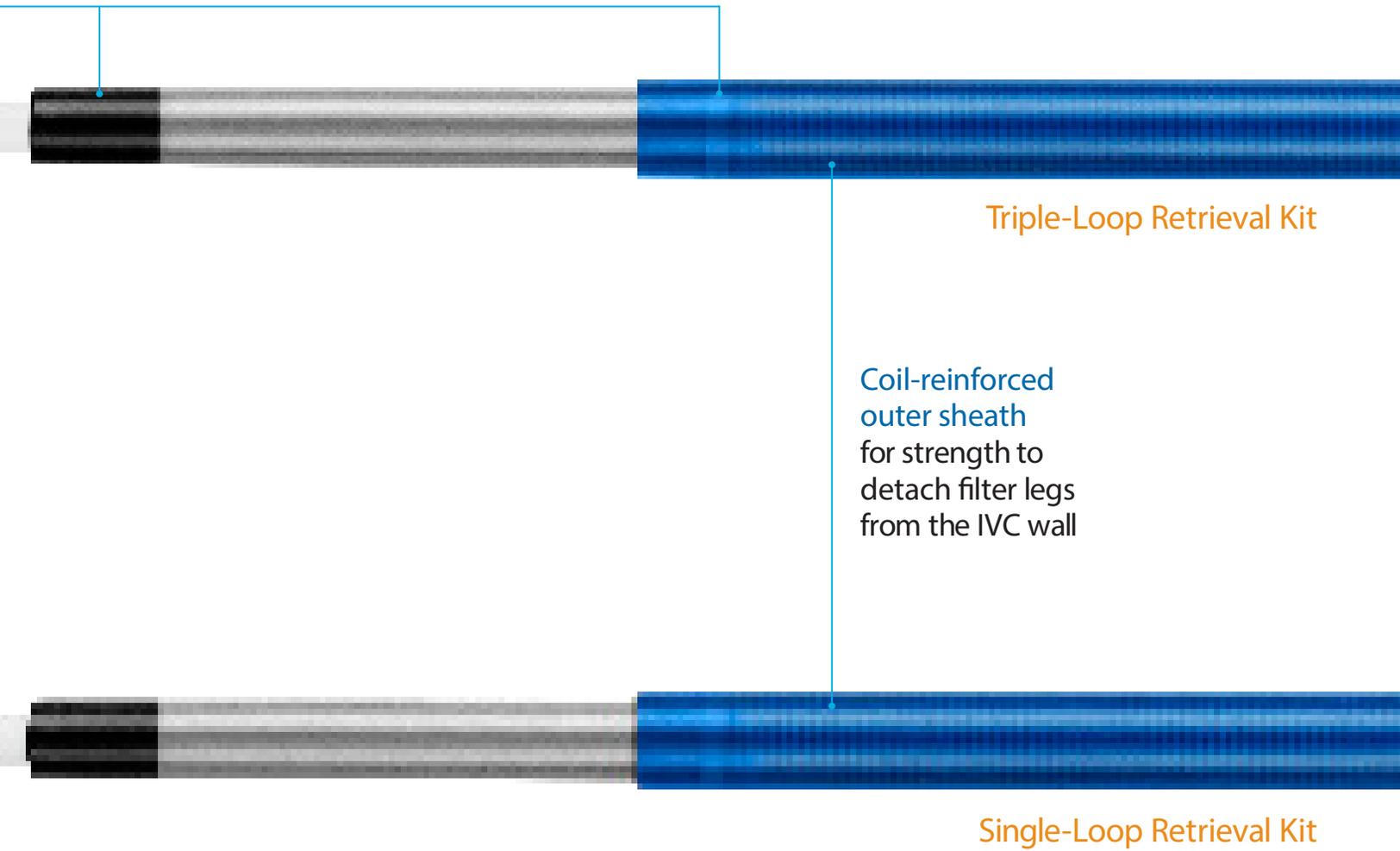

Triple-Loop Retrieval Kit

Coil-reinforced
outer sheath
for strength to
detach filter legs
from the IVC wall

Single-Loop Retrieval Kit


Ordering Information

| Single-Loop Retrieval Kit | | | |
|---------------------------|----------------|-----------------|----------------|
| CATALOG # | INNER SHEATH | OUTER SHEATH | SNARE CATHETER |
| 380911120 | 9F (ID) x 64cm | 11F (ID) x 60cm | 7F (OD) x 76cm |

| Triple-Loop Retrieval Kit | | | |
|---------------------------|----------------|-----------------|----------------|
| CATALOG # | INNER SHEATH | OUTER SHEATH | SNARE CATHETER |
| 380911330 | 9F (ID) x 64cm | 11F (ID) x 60cm | 7F (OD) x 76cm |

Each Kit Includes

- A. 11F (ID) x 60cm Outer Sheath
- B. 9F (ID) x 64cm Inner Sheath
- C. 8F (OD) x 75cm Dilator
- D. Hemostasis Valve with Sideport
- E. 3-Way High Pressure Stopcock
- F. 7F (OD) x 76cm Snare Catheter with Tuohy-Borst Y-Port Adapter, 30mm x 93cm Triple-Loop Snare with Torque Handle
- OR
- G. 7F (OD) x 76cm Snare Catheter with Tuohy-Borst Y-Port Adapter, 20mm x 93cm Single-Loop Snare with Torque Handle

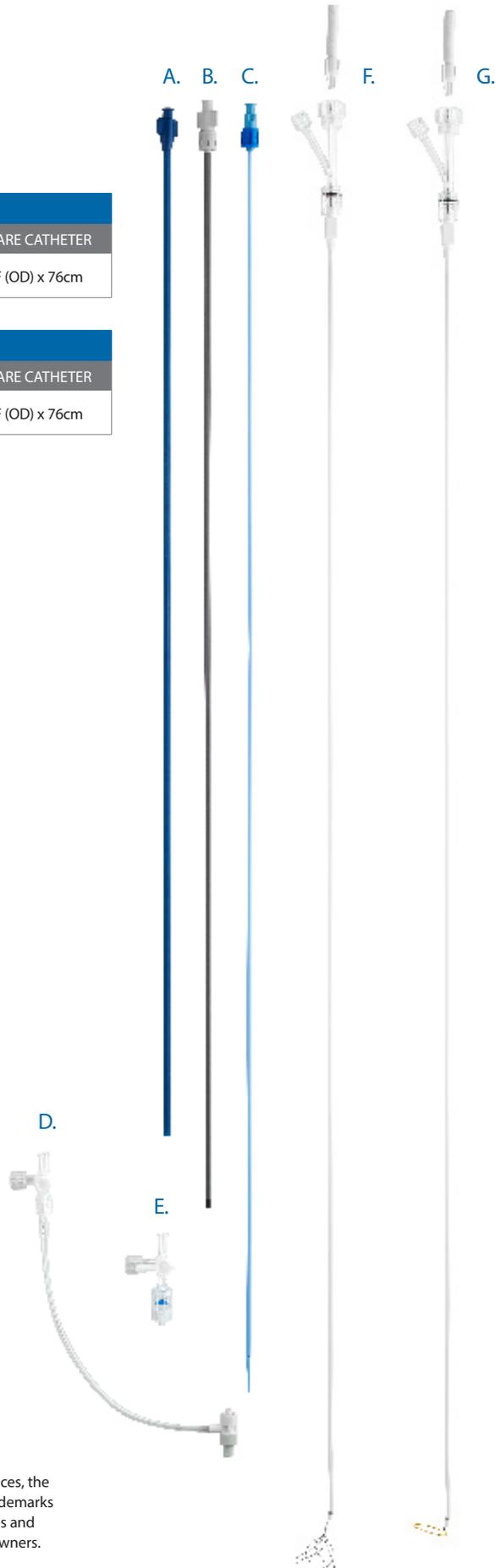

1. K070514, K152560, K073374, K112185

© 2019 Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries. All other trademarks are the property of their respective owners.

2600 Dallas Parkway, Suite 440, Frisco, TX 75034
 P: 800.927.4669 F: 903.677.9396 argonmedical.com

ARGON
 MEDICAL DEVICES

96-6500-09/0220B