

Tackle challenging retrievals with the CloverSnare™ 4-Loop Vascular Retriever



CloverSnare™ 4-LOOP VASCULAR RETRIEVER

High visibility

- The 6 Fr nitinol snare has a tantalum core for visibility.
- The radiopaque bands in each of the coaxial system's sheaths identify the distal tips.


Precise control

- Unique basket design maximizes torque control for retrieval.
- Slight angulation at catheter tip aids directional control to target.

More options

- Now you can retrieve filters, wires, coils, balloons, and more.



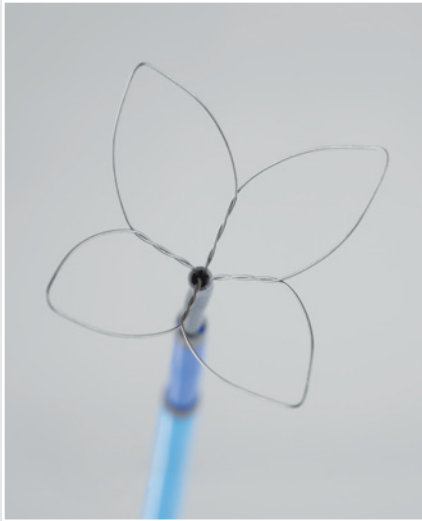

CloverSnare™

4-LOOP VASCULAR RETRIEVER

The CloverSnare 4-Loop Vascular Retriever is intended for use in the cardiovascular system to manipulate and retrieve foreign objects, including, but not limited to, wire guides, coils, balloons, catheters, and filters.

- The 6 Fr nitinol snare with tantalum core facilitates visibility and torque control.
- The radiopaque bands mark the distal ends of the 8 Fr inner and 10 Fr outer introducer sheaths.
- The hemostasis valve that is included on the inner sheath regulates flow.
- The introducer dilator, inner and outer coaxial introducer sheaths, and snare catheter lock together to provide a smooth transition.

Order Number	Reference Part Number	Snare Catheter Fr	Snare Catheter Length cm	Inner Sheath Length cm	Outer Sheath Length cm
G53008	VRS-6.0-90	6	90	85	80

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
 E-mail: cau.custserv@cookmedical.com



www.cookmedical.com