


# *Pinnacle<sup>®</sup> Destination<sup>®</sup>*

Guiding Sheath

GUIDE YOUR  
PROCEDURES  
THROUGH  
**CONSISTENT  
RELIABILITY**

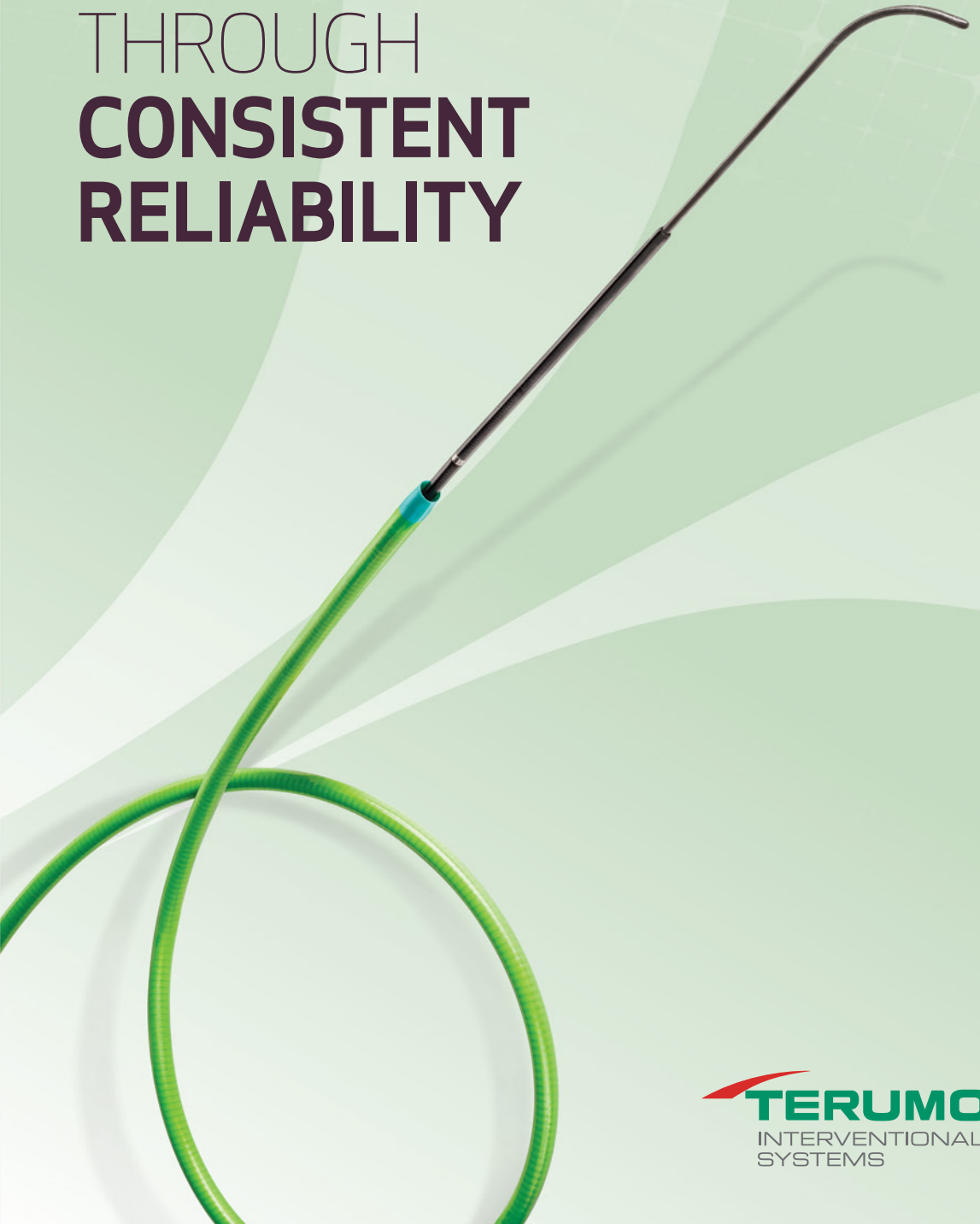

 **TERUMO**  
INTERVENTIONAL  
SYSTEMS



# PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you apply effective solutions and achieve **better outcomes for more patients.**

**TERUMO**  
INTERVENTIONAL  
SYSTEMS

# Pinnacle<sup>®</sup> Destination<sup>®</sup>

Guiding Sheath

## ALLOWS SMOOTH TRANSITION AND EASY PENETRATION

### Unique atraumatic tip<sup>1</sup>

- Minimizes potential for vessel damage

### Smooth transitions<sup>1</sup>

- Guidewire-to-dilator and dilator-to-sheath

### Provides lesion access<sup>1</sup>

- Exceptional trackability, even in the most challenging anatomy

## SHAPES AND SIZES TO MEET YOUR DIAGNOSTIC CHALLENGES

### Maintains shape and positioning<sup>1</sup>

- Outstanding durability, even in the longest procedures

### Enhances visualization for precise positioning<sup>1</sup>


- Triple radiopacity (sheath, dilator, and gold coil marker)

### Resists kinking<sup>1</sup>

- Stainless steel coil for consistent reliability

### Provides smooth movement<sup>1</sup>

- Coil reinforced tubing and a PTFE inner layer that minimizes friction

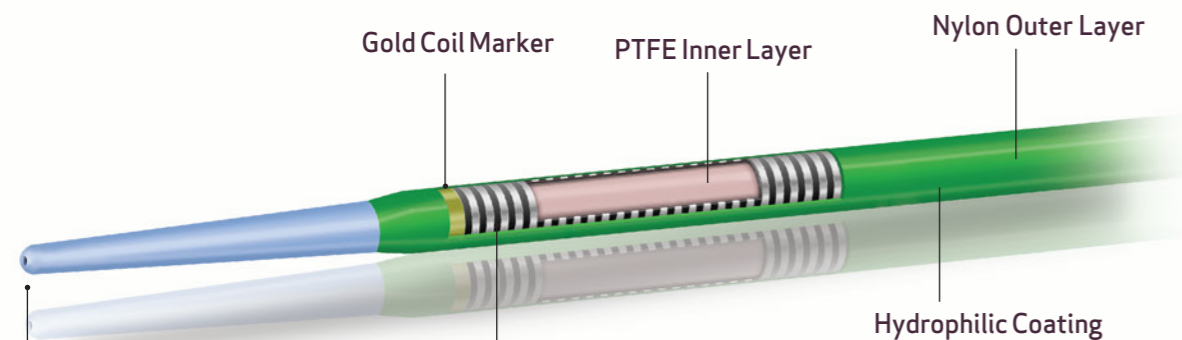

Stainless Steel Coil

Gold Coil Marker

PTFE Inner Layer

Nylon Outer Layer

Hydrophilic Coating

Detachable valve options

Nylon outer layer reinforced with stainless steel coil

Hydrophilic coating on distal end

Atraumatic tip

Cross Cut™ valve

Tuohy-Borst valve

### Accommodates anatomical variations in vessel access

- Straight
- Hockey Stick
- Multi-Purpose
- RDC (Renal Double Curve)
- LIMA (Left Internal Mammary Artery)


Designed to work with *Navicross<sup>®</sup>* and *Glidewire Advantage<sup>®</sup>*

Support Catheters

Peripheral Guidewires

# Pinnacle® Destination®

## Guiding Sheath

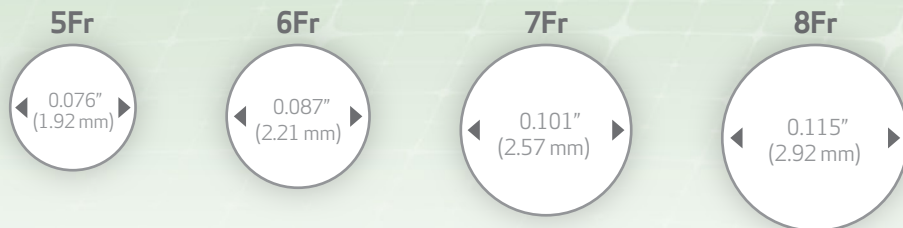

## PINNACLE® DESTINATION® Guiding Sheaths

Packed 1 per box and includes dilator. All dilators are 0.038" wire compatible.

French Size	Product Code	Length (cm)	Hydrophilic Coating Length (cm)	Valve Type	Curve Style	Curve Shape
5 Fr	54-54501	45	35	Cross-Cut	Straight	
	54-54502	45	5	Cross-Cut	Hockey Stick	
	54-54503	45	5	Cross-Cut	Multi-Purpose	
	54-54504	45	5	Cross-Cut	RDC	
6 Fr	54-64501	45	35	Cross-Cut	Straight	
	54-64506	45	35	Tuohy-Borst	Straight	
	RSR01	45	5	Cross-Cut	Straight	
	RSR02	45	5	Cross-Cut	Hockey Stick	
	RSR03	45	5	Cross-Cut	Multi-Purpose	
	RSR07	45	5	Tuohy-Borst	Straight	
	RSR08	45	5	Tuohy-Borst	Hockey Stick	
	RSR09	45	5	Tuohy-Borst	Multi-Purpose	
	RSR13	45	5	Cross-Cut	RDC	
	RSR14	45	5	Cross-Cut	LIMA	
	RSR17	45	5	Tuohy-Borst	RDC	
	RSR18	45	5	Tuohy-Borst	LIMA	
	RSP01	65	35	Cross-Cut	Straight	
	RSP03	65	35	Tuohy-Borst	Straight	
	RSC01	90	15	Tuohy-Borst	Straight	
	RSC03	90	15	Tuohy-Borst	Multi-Purpose	
RSC05	90	15	Cross-Cut	Straight		
RSC07	90	15	Cross-Cut	Multi-Purpose		
7 Fr	54-74501	45	35	Cross-Cut	Straight	
	54-74506	45	35	Tuohy-Borst	Straight	
	RSR04	45	5	Cross-Cut	Straight	
	RSR05	45	5	Cross-Cut	Hockey Stick	
	RSR06	45	5	Cross-Cut	Multi-Purpose	
	RSR10	45	5	Tuohy-Borst	Straight	
	RSR11	45	5	Tuohy-Borst	Hockey Stick	
	RSR12	45	5	Tuohy-Borst	Multi-Purpose	
	RSR15	45	5	Cross-Cut	RDC	
	RSR16	45	5	Cross-Cut	LIMA	
	RSR19	45	5	Tuohy-Borst	RDC	
	RSR20	45	5	Tuohy-Borst	LIMA	
	RSP02	65	35	Cross-Cut	Straight	
	RSP04	65	35	Tuohy-Borst	Straight	
	RSC02	90	15	Tuohy-Borst	Straight	
	RSC04	90	15	Tuohy-Borst	Multi-Purpose	
RSC06	90	15	Cross-Cut	Straight		
RSC08	90	15	Cross-Cut	Multi-Purpose		
8 Fr	54-84501	45	35	Cross-Cut	Straight	
	54-84506	45	35	Tuohy-Borst	Straight	
	54-86501	65	35	Cross-Cut	Straight	
	54-86506	65	35	Tuohy-Borst	Straight	
	54-89001	90	60	Cross-Cut	Straight	
	54-89006	90	60	Tuohy-Borst	Straight	

**FIND OUT MORE**



Phone: 800.862.4143



terumo.com



Fax: 800.411.5870

**RX ONLY.** Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

**Reference:**

1. Data on File.

©2019 Terumo Medical Corporation. All rights reserved.

All brand names are trademarks or registered trademarks of Terumo. PM-02093

