


OptiSeal® Valved Peelable Introducer

Smooth Introduction, Reliable Peel



Biphasic (two seal) valve designed to reduce risk of air embolism¹

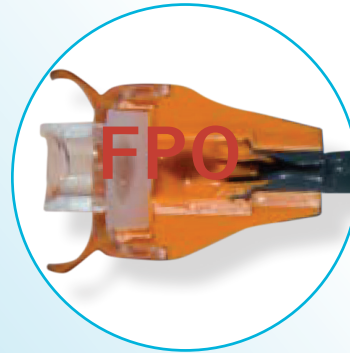


Symmetric, ergonomically designed wings allow for a smooth, unidirectional crack and peel

Optimal insertion and sealing during lead placement


Advanced safety and simplicity

- Designed to limit the clinician's exposure to blood-borne pathogens¹
- Smooth sheath-to-dilator transition facilitates atraumatic vessel access
- Color coded by French size for easy identification

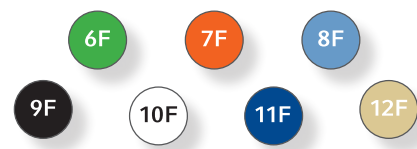

Kit contents

- 13cm or 25cm sheath & corresponding dilator
- 0.035", 60cm guidewire
- One 18G needle
- 12cc syringe

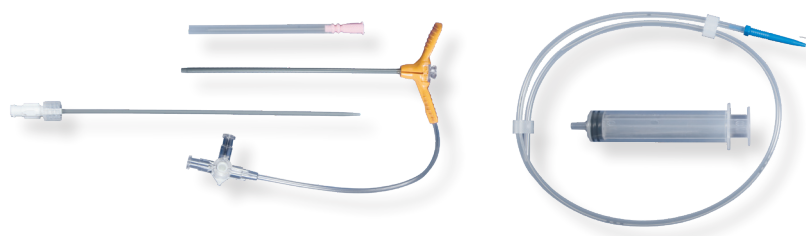

Ordering Information

OPTISEAL® VALVED PEELABLE INTRODUCER

MODEL	SIZE	SHEATH I.D.	SHEATH LENGTH	COLOR	BOX QUANTITY
1003975-001	6F	0.086"	13 cm	Dark Green	5
1003975-002	7F	0.099"	13 cm	Orange	5
1003975-003	8F	0.112"	13 cm	Light Blue	5
1003975-004	9F	0.125"	13 cm	Black	5
1003975-005	10F	0.138"	13 cm	White	5
1003975-006	11F	0.151"	13 cm	Blue	5
1003975-007	12F	0.164"	13 cm	Beige	5
1003975-008	6F	0.086"	25 cm	Dark Green	5
1003975-009	7F	0.099"	25 cm	Orange	5
1003975-010	8F	0.112"	25 cm	Light Blue	5
1003975-011	9F	0.125"	25 cm	Black	5
1003975-012	10F	0.138"	25 cm	White	5
1003975-013	11F	0.151"	25 cm	Blue	5
1003975-014	12F	0.164"	25 cm	Beige	5

*Individual units are not sold separately.

The OptiSeal Valved PTFE Peelable Introducer is intended for use in percutaneous insertion of pacing leads or catheters in the venous system.

References:

1. Vesely TM, Fazzaro AG, Gheradini D, Preliminary Evaluation of Valved Introducer Sheath for the Insertion of Tunneled Hemodialysis Catheters, Seminars in Dialysis, 2004; 17(1): 65-68

Please see instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

Teleflex, the Teleflex logo, and OptiSeal are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Other names may be trademarks of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 12/2020.

© 2020 Teleflex Incorporated. All rights reserved. MC-006089