

SPECTRALINE™ SYRINGES

A FULL SPECTRUM
OF CHOICES



Spectraline™ Syringes offer a choice of five different plunger colors and come packaged with a label set of eleven common injectate names for easy identification of fluid applications. They are available in 10ml and 20ml configurations for use in hand injecting saline, contrast, heparin or other fluids.

SPECTRALINE™ SYRINGES

- 10ml and 20ml syringes are available in five colors and come packaged with a set of eleven labels for common injectates such as nitro, local, lidocaine, saline, contrast, heparin and xylocaine.
- Clear polycarbonate barrel provides easy visualization of fluid levels and allows for accurate debubbling.
- The smooth operation of the piston through the polycarbonate barrel provides one handed operation and reduces hand fatigue.
- Large 0.100 inch I.D. syringe tip delivers injections with less force.
- Wide barrel flanges with finger grips prevent slippage and improve control.
- Ergonomic palm pad plunger provides optimum comfort for all size hands.
- Piston stop provides airtight connection and prevents piston from popping out of the barrel while aspirating.
- Double threaded male luer lock provides an airtight and secure connection.
- Rigid ABS piston maintains strength during forceful injections.

FOR MORE INFORMATION OR TO PLACE AN ORDER, CONTACT YOUR B. BRAUN INTERVENTIONAL SYSTEMS REPRESENTATIVE OR CALL 1-877-VENA CAV (836-2228)

Spectraline™ Syringes are indicated for use in cardiology and radiology procedures.


Refer to the Spectraline™ Syringes product insert for relevant warnings, precautions, complications and contraindications.

Rx only

CV827-Rev. 10/13

Spectraline™ is a trademark of B. Braun Medical Inc.

A FULL SPECTRUM OF CHOICES

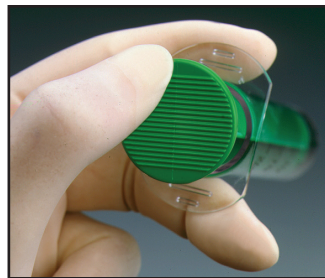

ORDERING INFORMATION: Spectraline™ Syringes

Plunger Color	Reorder Number	Syringe Volume	Units per Case
White	622600	10ml	25
Blue	622601	10ml	25
Green	622602	10ml	25
Yellow	622603	10ml	25
Red	622604	10ml	25
White	622610	20ml	25
Blue	622611	20ml	25
Green	622612	20ml	25
Yellow	622613	20ml	25
Red	622614	20ml	25

Each syringe comes individually packaged with a set of eleven labels.



Distributed by:
B. Braun Interventional Systems Inc.
 824 Twelfth Avenue
 Bethlehem, PA 18018 USA
 Tel: 1-877-VENA CAV (836-2228) (USA)
 Fax: 1-610-849-1334
 www.bisusa.org



Manufacturer:
Merit Medical Systems, Inc.
 1600 West Merit Parkway
 South Jordan, Utah 84095