

The Power to Restore Flow.

CLEANERXT™ and CLEANER15™ Rotational Thrombectomy Systems




ARGON
MEDICAL DEVICES

The Argon Medical Devices' family of thrombectomy products




Combining **Torque** with **Trackability**

- **NEW** – Dual power cells increase duration and run time, maintaining torque through tortuosity
- **NEW** – Improved drive shaft design increases efficiency of power transfer and torque
- 6F low profile design enables excellent steerability through vessel tortuosity

The **POWER** to Restore Flow


- Increased sinusoidal wire amplitude
– From 9 mm to 15 mm
- Increased wire diameter
– From .035” to .044”
- Dual power cells maintain torque through tortuosity
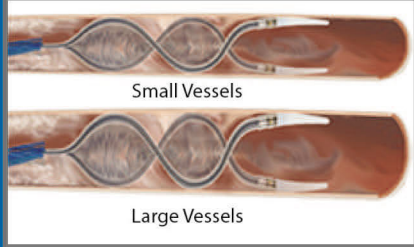
| Power Comparison | Power Cells | Voltage | RPM | Conforms to lumen diameter to clean wall-adherent thrombus  |
|-----------------------|-------------|---------|------|---|
| CLEANERXT | 2 | 9 Volts | 4000 | |
| CLEANER15 | 2 | 9 Volts | 4000 | |
| Arrow Tretotola™ PTD® | 1 | 6 Volts | 3000 | |

CLEANERXT and CLEANER15 are indicated for mechanical declotting of native vessel dialysis fistulae and synthetic dialysis access grafts


Innovative Sinusoidal Wire Design

- Radiopaque S-wire conforms to varying lumen diameters, actively cleaning wall-adherent thrombus
- Fluid vortex gently macerates thrombus while reducing the risk of endothelial damage¹
- S-wire is capable of cleaning around opposing sheaths reducing procedure time and sheath manipulation

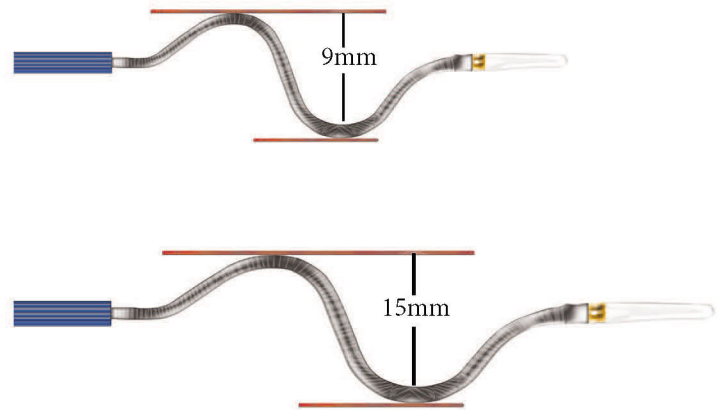

Atraumatic Mechanical Thrombectomy


- Safe and effective in both AV fistulae and grafts
- Guidewire-like design provides excellent steerability
- Available in 65 cm and 135 cm catheter lengths
- 3-way sideport and distal side hole on catheter allows for infusion of fluids and contrast media

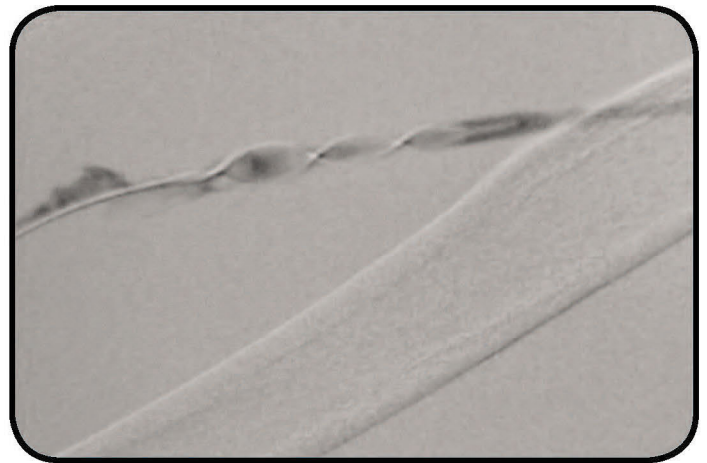

Used with permission of Mark Goldberg, MD

Ergonomic One Piece Design

- Procedure ready – no assembly required
- Intuitive to operate
- Quiet operation allows for patient comfort during procedure
- Completely disposable – no capital equipment required

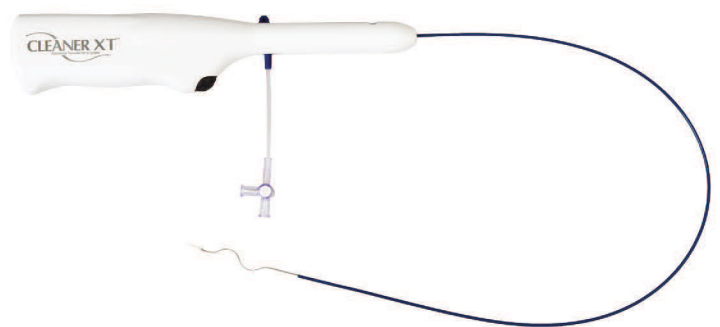

For more information visit www.argonmedical.com



CLEANER XT ORDERING INFORMATION

| QUANTITY | CATALOG NUMBER | CATHETER OD | LENGTH | SINUSOIDAL WIRE AMPLITUDE |
|----------|----------------|-------------|--------|---------------------------|
| _____ | 700009XT | 6F | 65cm | 9mm |
| _____ | 700109XT | 6F | 135cm | 9mm |



CLEANER 15 ORDERING INFORMATION

| QUANTITY | CATALOG NUMBER | CATHETER OD | LENGTH | SINUSOIDAL WIRE AMPLITUDE |
|----------|----------------|-------------|--------|---------------------------|
| _____ | 700015 | 7F | 65cm | 15mm |
| _____ | 700115 | 7F | 135cm | 15mm |

PHYSICIAN'S SIGNATURE



1. F. Vega, Safety Comparison of Cleaner Rotational Thrombectomy System versus Arrow PTD in Native Vessels

Argon Medical Devices, Inc. own, use or have applied for the following trademarks or service marks: CLEANERXT, CLEANER15, Argon Medical Devices, Inc. All other trademarks are trademarks of their respective owners or holder.

1445 Flat Creek Road, Athens, TX 75751
P: 800.927.4669 F: 903.677.9396 argonmedical.com

96-6000-37/0517B
DDM/SS
Copyright 2017 Argon Medical Devices, Inc.
Printed in USA