

DuraFlow 2


CHRONIC DIALYSIS CATHETER



DuraFlow 2

CHRONIC DIALYSIS CATHETER

Engineered for Durability

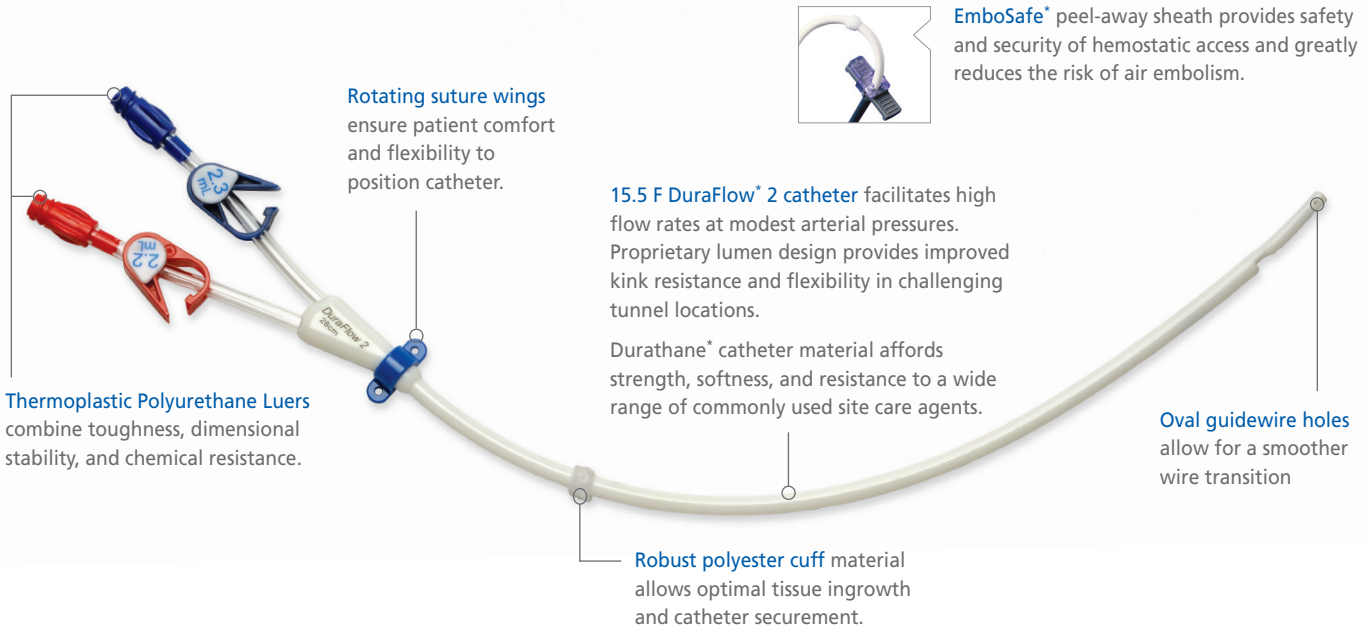

Rotating suture wings ensure patient comfort and flexibility to position catheter.

Thermoplastic Polyurethane Luers combine toughness, dimensional stability, and chemical resistance.

15.5 F DuraFlow² catheter facilitates high flow rates at modest arterial pressures. Proprietary lumen design provides improved kink resistance and flexibility in challenging tunnel locations.

Durathane[®] catheter material affords strength, softness, and resistance to a wide range of commonly used site care agents.

Oval guidewire holes allow for a smoother wire transition

Robust polyester cuff material allows optimal tissue ingrowth and catheter securement.

EmboSafe[®] peel-away sheath provides safety and security of hemostatic access and greatly reduces the risk of air embolism.

Ordering Information (5 units per box)

Description	Basic Kit	VascPak Kit	Catheter Only
20 cm	<input type="checkbox"/> H787103012145	—	—
22 cm	<input type="checkbox"/> H787103012125	—	—
24 cm	<input type="checkbox"/> H787103012015	<input type="checkbox"/> H787103012185	<input type="checkbox"/> H787103017015
28 cm	<input type="checkbox"/> H787103012025	<input type="checkbox"/> H787103012195	<input type="checkbox"/> H787103017025
32 cm	<input type="checkbox"/> H787103012035	<input type="checkbox"/> H787103012205	<input type="checkbox"/> H787103017035
36 cm	<input type="checkbox"/> H787103012045	<input type="checkbox"/> H787103012215	<input type="checkbox"/> H787103017045
40 cm	—	—	<input type="checkbox"/> H787103017055
48 cm	—	—	<input type="checkbox"/> H787103017155
55 cm	—	—	<input type="checkbox"/> H787103017145

Ordering Information (1 unit per box)

Description	Basic Kit	VascPak Kit
24 cm	<input type="checkbox"/> H787103012011	—
28 cm	<input type="checkbox"/> H787103012021	—
32 cm	<input type="checkbox"/> H787103012031	—
36 cm	<input type="checkbox"/> H787103012041	—
40 cm	<input type="checkbox"/> H787103012055	<input type="checkbox"/> H787103012225
48 cm	<input type="checkbox"/> H787103012155	<input type="checkbox"/> H787103012235
55 cm	<input type="checkbox"/> H787103012115	<input type="checkbox"/> H787103012245

BASIC KIT INCLUDES: (1) DuraFlow 2 Dialysis Catheter, (1) 16F Valved Breakaway Introducer, (1) 12F Dilator, (1) 14F Dilator, (1) Tri-Ball Tunneler with Sleeve, (2) Injection Caps, (1) #11 Blade Safety Scalpel, (1) 0.038" J/Flex Guidewire, (1) Adhesive Wound Dressing, (1) 18 Ga. x 2-3/4" Introducer Needle

VASCPAK[®] KIT INCLUDES: (1) DuraFlow 2 Dialysis Catheter, (1) Adhesive Wound Dressing, (1) 16F Valved Peelable Sheath, (1) Tri-ball Tunneler with Sleeve, (2) Injection Caps, (1) 12F Vessel Dilator, (1) 14F Vessel Dilator

IMPORTANT RISK INFORMATION

INDICATIONS FOR USE: The DuraFlow 2 Hemodialysis Catheter is indicated for use in attaining long-term vascular access for hemodialysis and apheresis. It may be inserted percutaneously and is primarily placed in the internal jugular vein of an adult patient. Alternate insertion sites include subclavian vein as required. Catheters greater than 40 cm are intended for femoral vein insertion. This catheter is indicated for > 30 days long term placement. The valved peelable introducer is

intended for use in the percutaneous insertion of catheters in the venous system.

CONTRAINDICATIONS: The catheter is intended for long-term vascular access only and should not be used for any purpose other than indicated in these instructions. The valved peelable introducer sheath is not designed for use in the arterial system.

WARNINGS AND PRECAUTIONS: Please see package insert for complete list of warnings and precautions

POTENTIAL COMPLICATIONS: Consult package insert for a complete list of potential complications.

CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician.

Indications, contraindications, warnings and instructions for use can be found in the instructions for use supplied with each device. Observe all instructions prior to use. Failure to do so may result in patient complications.



USA > 14 Plaza Drive, Latham, NY 12110 > tel: 800-772-6446 or 518-798-1215 > fax: 518-798-1360

International > Haaksbergweg 75 (Margrietoren), 1101 BR, Amsterdam Z-O > The Netherlands

tel: +31 (0)20 753 2949 > fax: +31 (0)20 753 2939

www.angiodynamics.com

*AngioDynamics, the AngioDynamics logo, DuraFlow, EmboSafe, Durathane and Vascpak are trademarks and/or registered trademarks of AngioDynamics, Inc., an affiliate or a subsidiary.

© 2015 AngioDynamics, Inc. MLC 456 US Rev E 04/15