

Angioplasty Inflation Device

Designed to take the pressure



Provides precise control with repeatable, accurate pressure for use in the most demanding interventional procedures.

Angioplasty Inflation Device

Accurate. Responsive. Consistent.

Highlights:

- 30 ATM luminescent analog pressure gauge accommodates most procedures and registers inflation pressure immediately without decay or drift.
- 25cc clear polycarbonate syringe barrel provides rapid deflation, accommodates large balloon sizes and permits easy debubbling.
- O-ring syringe piston design provides accurate pressure and minimizes pressure decay.
- Winged locking mechanism permits easy deflation at high pressures and engages more threads for stability under high pressure.


For more information or to place an order, contact your B. Braun Interventional Systems Inc. representative or call 1-877-VENA CAV (836-2228)

Ordering Information:

Inflation Device		
Product Description	REF Number	Case Qty.
Inflation Device	622510	10
Inflation Device Kit with single Y-connector and plastic insertion tool	622511	10

Single (9.5 F lumen) Y-connector		
Single Y-connector with metal insertion tool	610400	5
Single Y-connector with plastic insertion tool	610402	5

Y-connector Kits (9.5 F lumen)		
Single Y-connector, metal insertion tool and torque device	610420	5
Single Y-connector, plastic insertion tool and torque device	610403	5

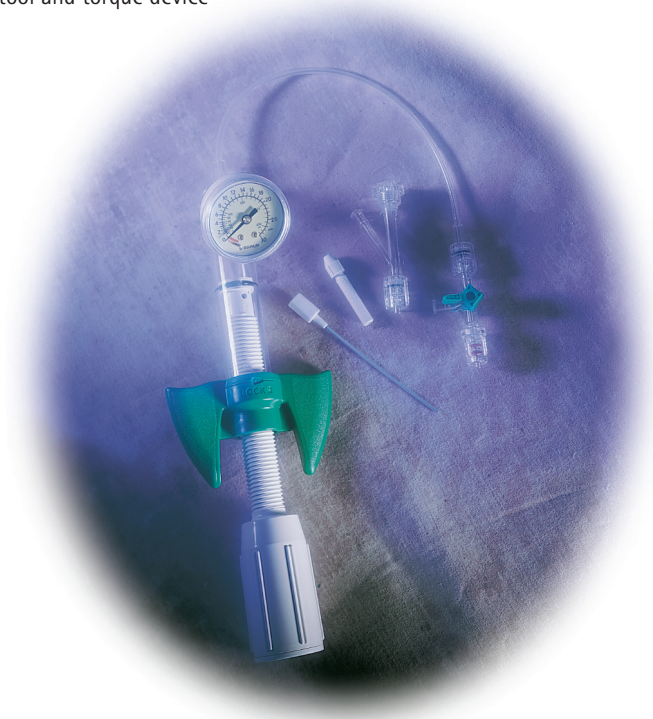

Manufactured by:
B. Braun Interventional Systems Inc.
 824 Twelfth Avenue
 Bethlehem, PA 18018 USA
 Tel: 1-877-VENA CAV (836-2228)
 Fax: 610-266-3982
www.bbraunusa.com