


LeMaitre® Embolectomy Catheters

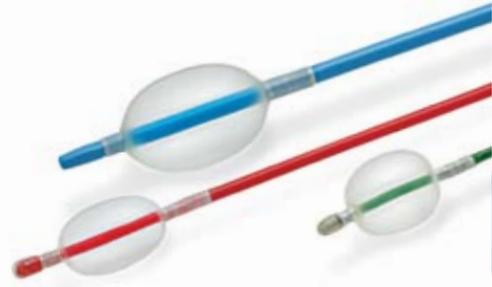

LeMaitre® Embolectomy Catheters

ORDERING INFORMATION

Ballon Diameter	Ballon Volume	Diameter	Color	Length	Guidewire	Model #
Over-the-Wire Embolectomy Catheter Compatible with standard guidewire sizes. Stopcock and syringe included.						
6 mm	0.20 ml	3F	Green	40 cm	.018"	1651-34
6 mm	0.20 ml	3F	Green	80 cm	.018"	1651-38
10 mm	0.75 ml	4F	Red	40 cm	.025"	1651-44
10 mm	0.75 ml	4F	Red	80 cm	.025"	1651-48
12 mm	1.50 ml	5F	White	40 cm	.028"	1651-54
12 mm	1.50 ml	5F	White	80 cm	.028"	1651-58
12 mm	1.50 ml	5F Plus	White	40 cm	.035"	1651-84
12 mm	1.50 ml	5F Plus	White	80 cm	.035"	1651-88
13 mm	1.60 ml	6F	Blue	40 cm	.035"	1651-64
13 mm	1.60 ml	6F	Blue	80 cm	.035"	1651-68
14 mm	1.75 ml	7F	Yellow	80 cm	.038"	1651-78

Single Lumen Embolectomy Catheter

4.5 mm	0.05 ml	2F	Purple	40 cm		1601-24
4.5 mm	0.05 ml	2F	Purple	60 cm		1601-26
8.0 mm	0.20 ml	3F	Green	40 cm		1601-34
8.0 mm	0.20 ml	3F	Green	80 cm		1601-38
10.5 mm	0.75 ml	4F	Red	40 cm		1601-44
10.5 mm	0.75 ml	4F	Red	80 cm		1601-48
13.0 mm	1.50 ml	5F	White	80 cm		1601-58
13.5 mm	1.60 ml	6F	Blue	80 cm		1601-68
14.0 mm	1.75 ml	7F	Yellow	80 cm		1601-78

NovaSil® Latex-Free Single Lumen Embolectomy Catheter

4 mm	0.05 ml	2F	Purple	40 cm		1801-24
4 mm	0.05 ml	2F	Purple	60 cm		1801-26
6 mm	0.20 ml	3F	Green	40 cm		1801-34
6 mm	0.20 ml	3F	Green	80 cm		1801-38
9 mm	0.60 ml	4F	Red	40 cm		1801-44
9 mm	0.60 ml	4F	Red	80 cm		1801-48
11 mm	1.00 ml	5F	White	80 cm		1801-58
13 mm	1.60 ml	6F	Blue	80 cm		1801-68
14 mm	1.75 ml	7F	Yellow	80 cm		1801-78

LeMaitre Catheter Rack

Catheter Rack	A3000-00
Wall Mount	A3000-01

These specifications are not intended as a warranty. In the interest of product improvement, these specifications may be changed from time to time without notice. Please consult your sales representative for details.



Your Peripheral Vision™

LeMaitre and NovaSil are registered trademarks of LeMaitre Vascular, Inc.
©2009 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA. M0005 Rev. G 09/09

LeMaitre Vascular, Inc.

63 Second Avenue
Burlington, MA 01803
USA

T 781 221 2266

F 781 221 2223

LeMaitre Vascular GmbH

Otto-Volger-Str. 5a/b
65843 Sulzbach/Ts.
Germany

T +49 (0)6196 659230

F +49 (0)6196 527072

LeMaitre Vascular GK

1F Kubodera Twin Tower Bldg.
2-9-4 Kudan-minami, Chiyoda-ku
Tokyo 102-0074
Japan

T +81 (0)3 5215 5681

F +81 (0)3 5215 5682

5F Plus featuring two radiopaque bands for visibility



For More Information:

Contact your local LeMaitre
Vascular Sales Representative.

For Customer Service:

Call 800 628 9470

www.lemaitre.com



Your Peripheral Vision™



LeMaitre® Embolectomy Catheters

A COMPLETE RANGE OF QUALITY CATHETERS:

LeMaitre offers an expanded line of premium quality balloon catheters for rapid removal of emboli and thrombi. Available in Single Lumen, Over-the-Wire, and Latex-Free models, LeMaitre Catheters are:

Radiopaque.

For easy visibility of instrument under fluoroscopy

Pliable.

Eases introduction into difficult anatomy

Sterile.

Three-year sterility for latex balloons; five-year for latex-free

OVER-THE-WIRE

Guidance.


Compatible with standard guidewire to direct catheter to proper position

Dye Enhancement.

Add contrast media to facilitate emboli/thrombi imaging and improve fluoroscopy control

Irrigation.

Lubricates during catheter withdrawal

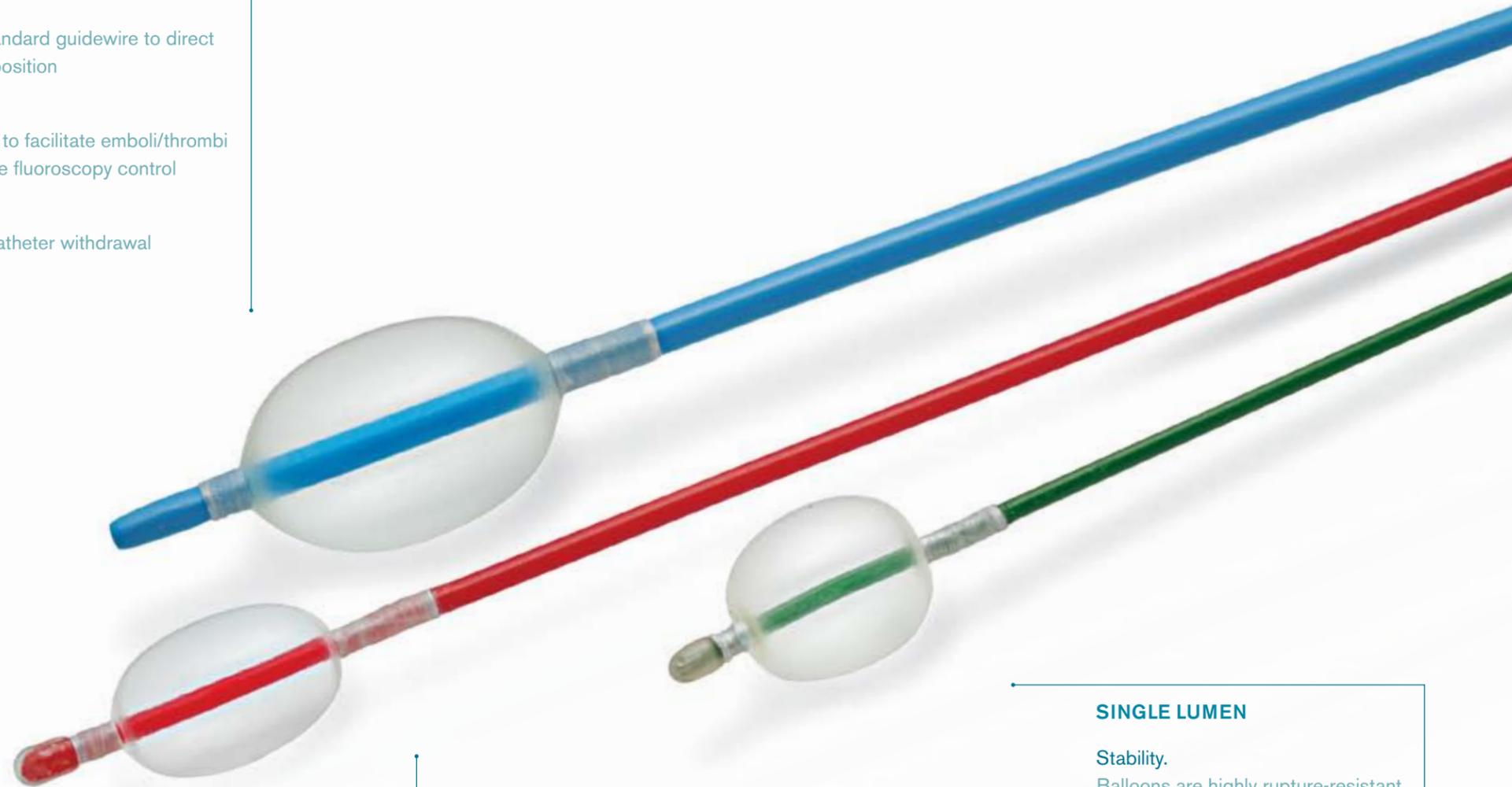

SINGLE LUMEN

Stability.

Balloons are highly rupture-resistant

Concentricity.

Balloons are well-centered

Progressive Inflation.

Balloons inflate at consistent rate as volume increases

LeMaitre also offers NovaSil® **LATEX-FREE** Catheters. This line of premium quality silicone balloon catheters provides a latex-free option for latex sensitivity.

5F PLUS



- Two radiopaque bands at proximal and distal ends of balloon for enhanced visibility under fluoroscopy
- Compatible with .035" guidewire and 6F introducer