



Navigate Complex  
Interventions at Every Turn

# Keep Driving Forward

The portfolio of Turnpike® Catheters places the power in your hands. Navigate tortuous anatomy with agility and precision. No matter the challenge, Turnpike Catheters are designed to help simplify the journey.



*"Turnpike® Catheters have a vast utility. They've given me the confidence to tackle a wide variety of highly complex lesions."\**


**Jon Robken, MD**  
Cardiovascular Medicine, P.C.

## Ideal Combination of Flexibility and Torque Response

Turnpike® Catheters feature a unique five-layer shaft construction that provides flexibility, torque and tracking over a 0.014" guidewire in complex coronary interventions.

## Outmaneuver Roadblocks

- + Complex lesions
- + Extreme tortuosity
- + Guidewire support

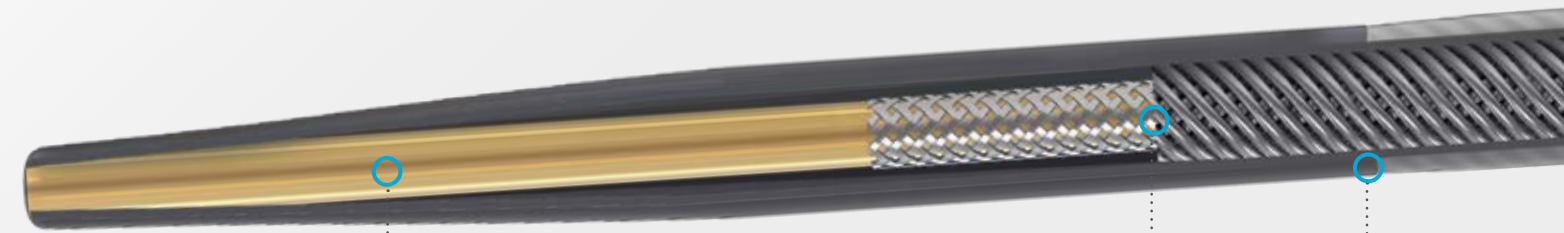

### Enhanced Guidewire Movement

An inner PTFE liner that spans the entire catheter to ensure smooth guidewire interaction throughout the case

### Responsive and Kink Resistant

Dual-layer coils and a braid make for an ideal combination of flexibility and torque response

### Exceptional Delivery

A polymer outer layer, paired with a hydrophilic coating on the distal 60 cm, provides excellent lubricity for delivery

### Turnpike® LP

*Dual coil tapers down to single coil 21 cm from distal tip for lower profile and increased flexibility in the distal shaft and tip*

\* This statement reflects the personal experience and opinion of the physician.





### Standard version, flexible tip

Workhorse antegrade device, as well as retrograde escalation device for increasing torque power and guidewire support<sup>1</sup> in resistant septal collaterals



5-layer shaft construction includes dual-layer coils and a braid for an ideal combination of flexibility and torque response

### Superior Bidirectional Torque Response\*

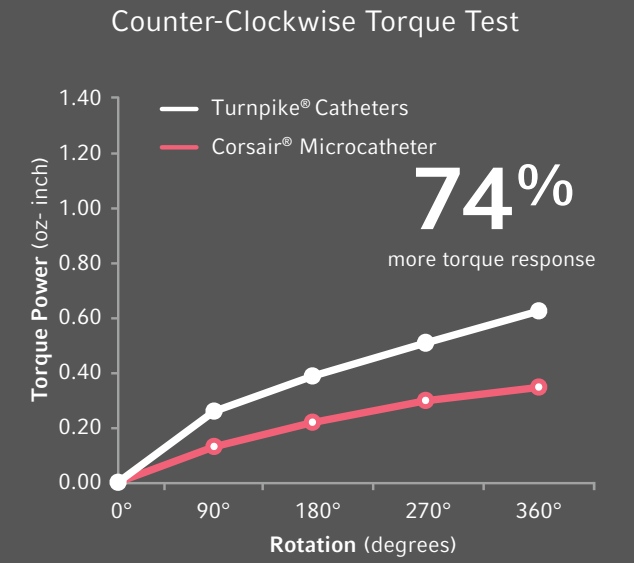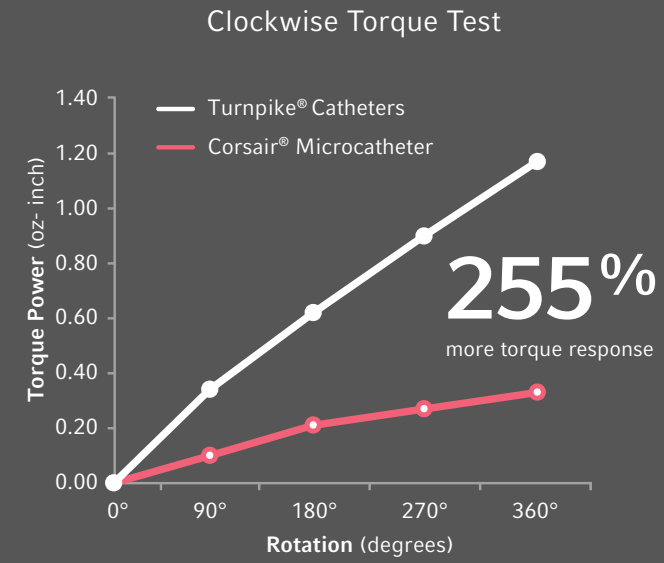


\* All values based on bench test data averages, n=3, performed by Teleflex. Bench test results may not necessarily be indicative of clinical performance. Data on file.

### Ordering Information

MODEL	VERSION	WORKING LENGTH	PROXIMAL SHAFT O.D.	DISTAL SHAFT O.D.	DISTAL TIP O.D.	COMPATIBLE GUIDEWIRE DIAMETER	COMPATIBLE GUIDE CATHETER SIZE
5642	Turnpike	135 cm	2.9 Fr. 0.97 mm / 0.038"	2.6 Fr. 0.86 mm / 0.034"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.
5643	Turnpike	150 cm	2.9 Fr. 0.97 mm / 0.038"	2.6 Fr. 0.86 mm / 0.034"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.

Packaged in quantities of 1 unit per box.

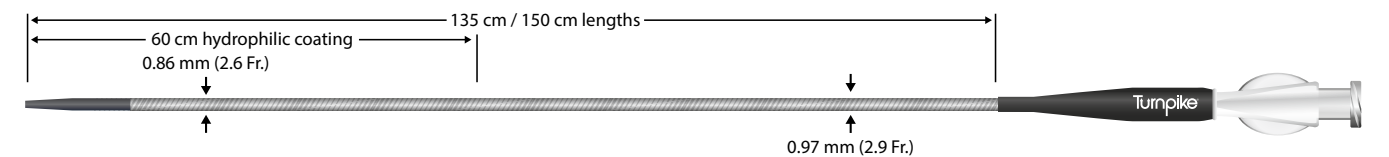


<sup>1</sup> Compared to Turnpike® LP Catheter. Data on file.

# Turnpike<sup>®</sup> Spiral catheter

## Spiral shaft, flexible tip

Front-line device for complex antegrade interventions

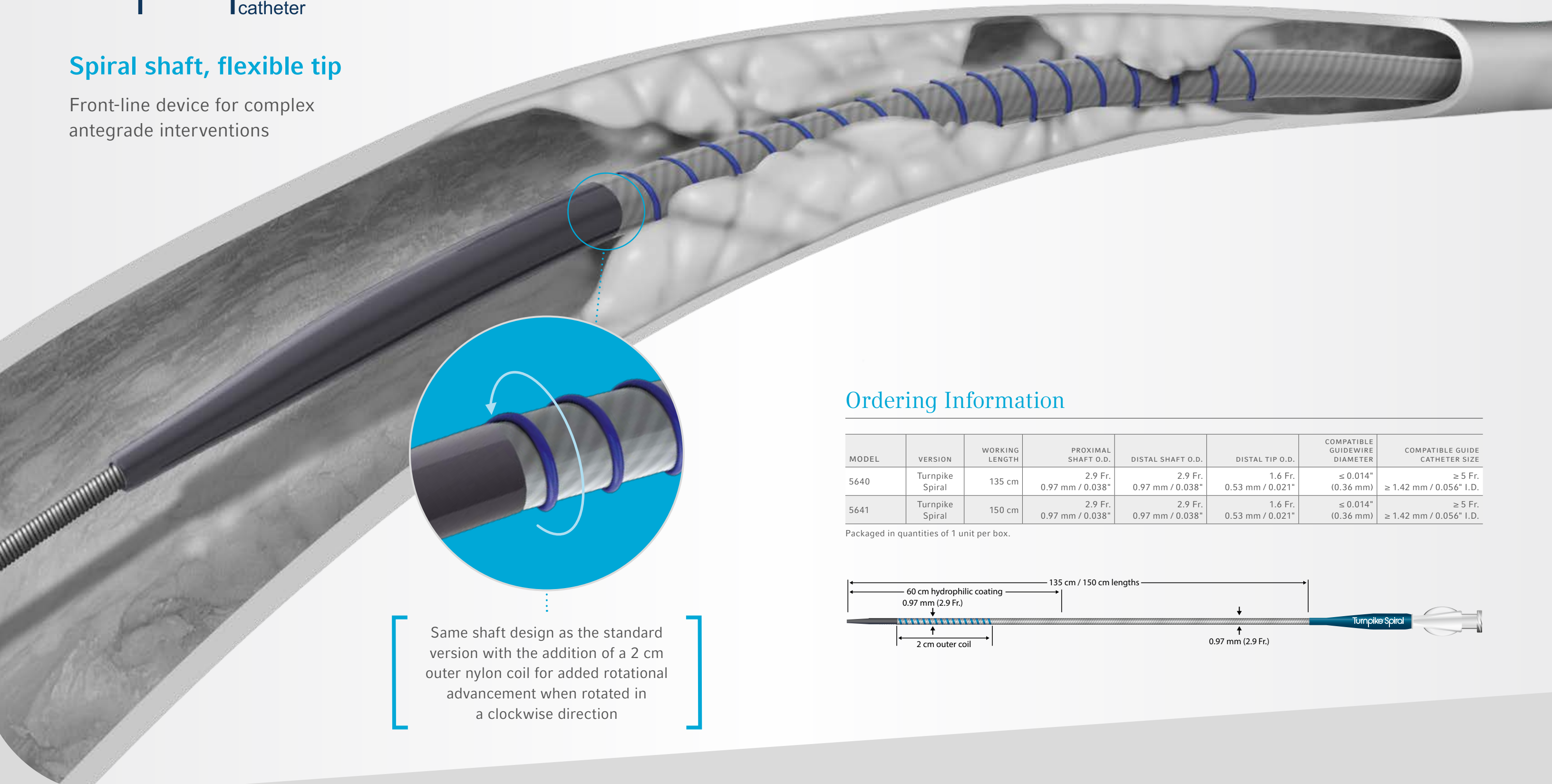


Same shaft design as the standard version with the addition of a 2 cm outer nylon coil for added rotational advancement when rotated in a clockwise direction

## Ordering Information

MODEL	VERSION	WORKING LENGTH	PROXIMAL SHAFT O.D.	DISTAL SHAFT O.D.	DISTAL TIP O.D.	COMPATIBLE GUIDEWIRE DIAMETER	COMPATIBLE GUIDE CATHETER SIZE
5640	Turnpike Spiral	135 cm	2.9 Fr. 0.97 mm / 0.038"	2.9 Fr. 0.97 mm / 0.038"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.
5641	Turnpike Spiral	150 cm	2.9 Fr. 0.97 mm / 0.038"	2.9 Fr. 0.97 mm / 0.038"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.



Packaged in quantities of 1 unit per box.

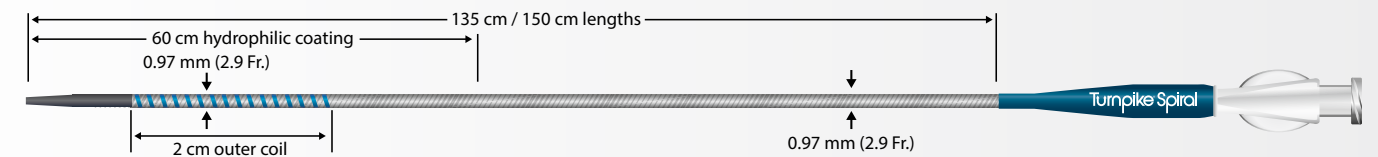

# Turnpike<sup>®</sup> Gold catheter

## Spiral shaft, threaded metal tip

Antegrade escalation device for resistant lesions

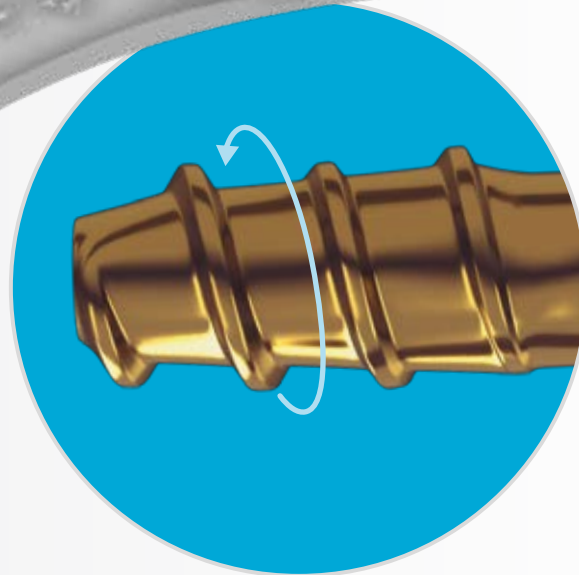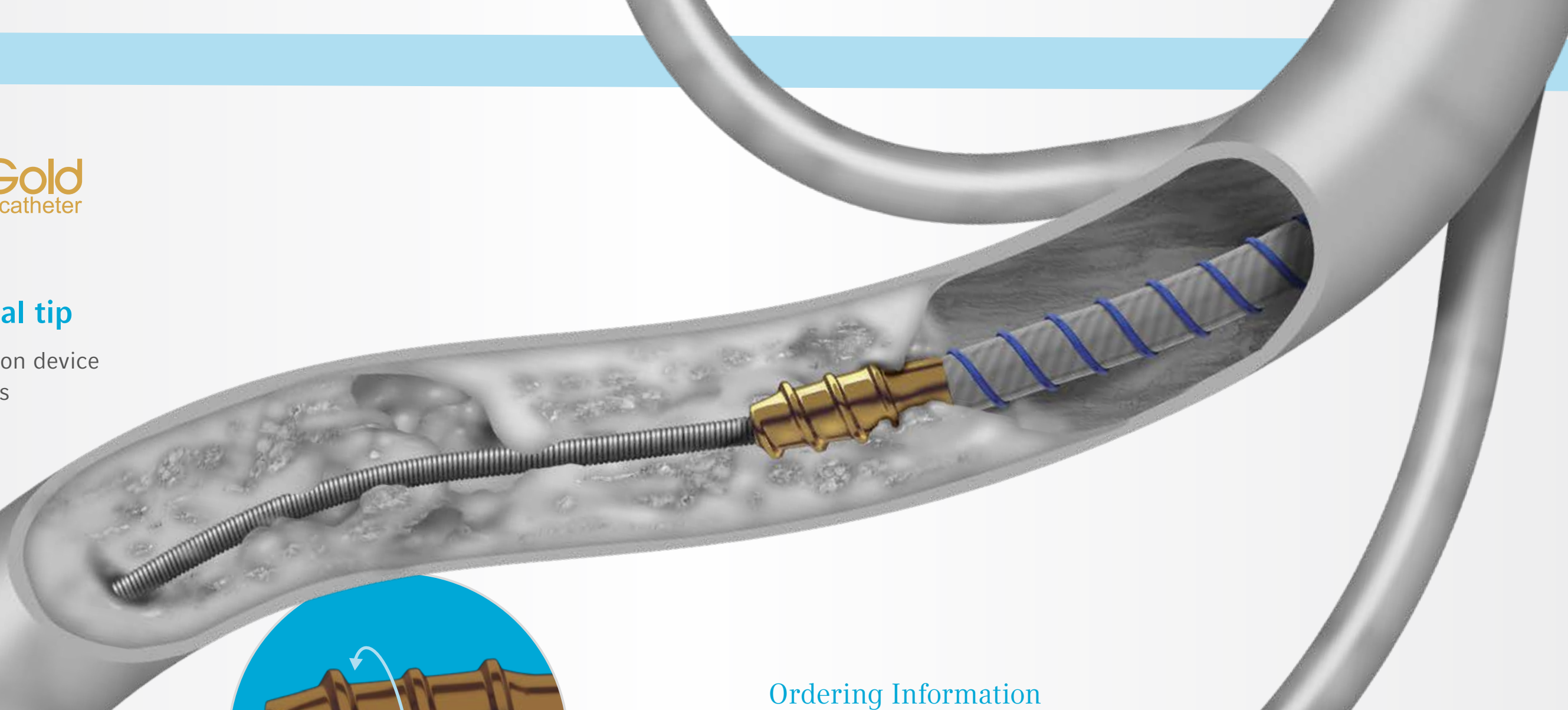


Gold-plated threads on distal tip provide leading edge rotational advancement with visibility under fluoroscopy

## Ordering Information

MODEL	VERSION	WORKING LENGTH	PROXIMAL SHAFT O.D.	DISTAL SHAFT O.D.	DISTAL TIP O.D.	COMPATIBLE GUIDEWIRE DIAMETER	COMPATIBLE GUIDE CATHETER SIZE
5621	Turnpike Gold	135 cm	2.9 Fr. 0.97 mm / 0.038"	2.9 Fr. 0.97 mm / 0.038"	2.1 Fr. 0.71 mm / 0.028"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.
5622	Turnpike Gold	150 cm	2.9 Fr. 0.97 mm / 0.038"	2.9 Fr. 0.97 mm / 0.038"	2.1 Fr. 0.71 mm / 0.028"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.


Packaged in quantities of 1 unit per box.

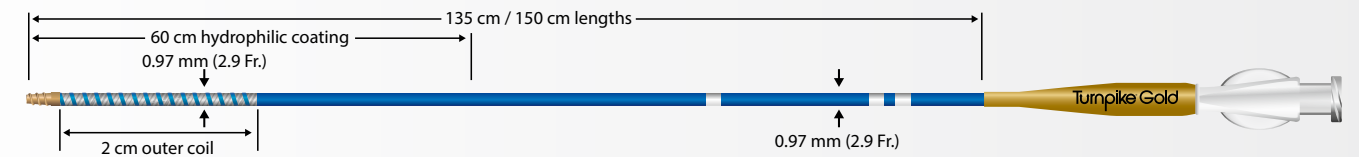

# Turnpike<sup>®</sup> LP catheter

## Low-profile shaft, flexible tip

Front-line device for retrograde interventions

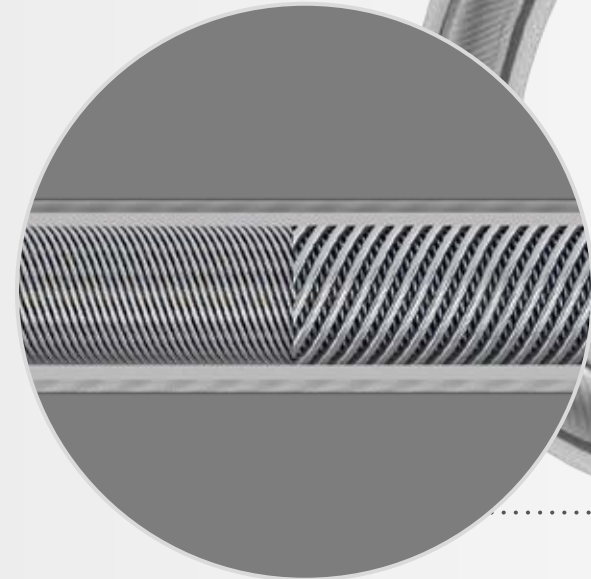


Distal Shaft Flexibility  
**79.4%**  
more flexible<sup>1</sup>


Tip Flexibility  
**85.7%**  
more flexible<sup>1</sup>

Dual-to-single coil transition 21 cm from distal tip for lower profile and increased flexibility in the distal shaft and tip


Turnpike<sup>®</sup> LP Catheter offers impressive torqueability with a smaller crossing profile than the Corsair<sup>®</sup> Microcatheter

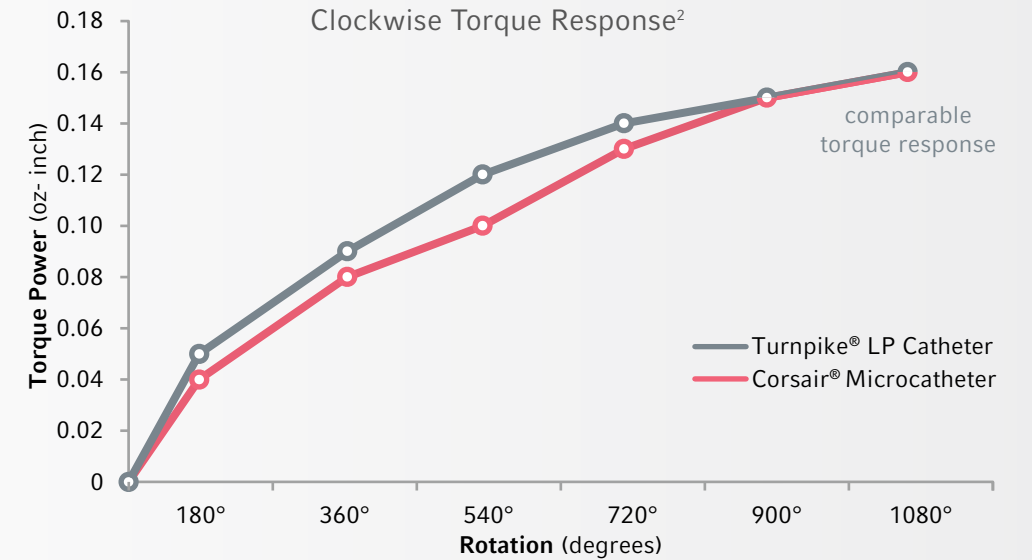

### Torqueable Microcatheters

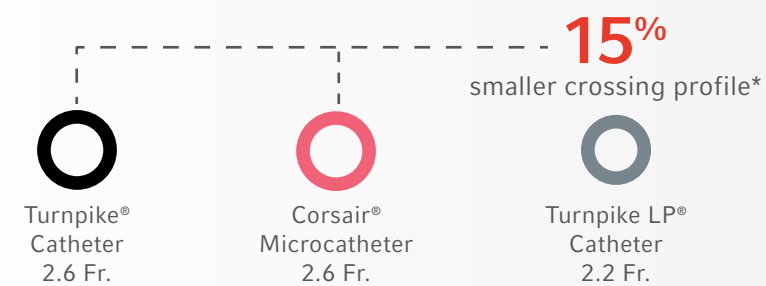

### Non-Torqueable Microcatheters

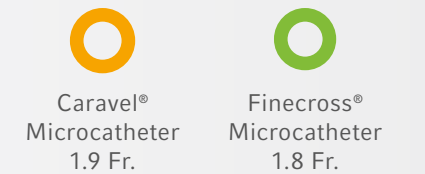


\* Based on diameter values reported on competitive websites at the time of this publication.

## Ordering Information

MODEL	VERSION	WORKING LENGTH	PROXIMAL SHAFT O.D.	DISTAL SHAFT O.D.	DISTAL TIP O.D.	COMPATIBLE GUIDEWIRE DIAMETER	COMPATIBLE GUIDE CATHETER SIZE
5638	Turnpike LP	135 cm	2.9 Fr. 0.97 mm / 0.038"	2.2 Fr. 0.74 mm / 0.029"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.
5639	Turnpike LP	150 cm	2.9 Fr. 0.97 mm / 0.038"	2.2 Fr. 0.74 mm / 0.029"	1.6 Fr. 0.53 mm / 0.021"	≤ 0.014" (0.36 mm)	≥ 5 Fr. ≥ 1.42 mm / 0.056" I.D.

Packaged in quantities of 1 unit per box.

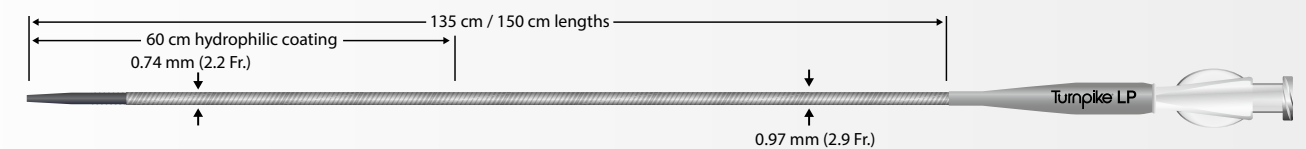

1 All values based on bench test data averages compared to the standard Turnpike Catheter, n=15, performed by Teleflex. Bench test results may not necessarily be indicative of clinical performance. Data on file.  
2 All values based on bench test data averages, n=2, performed by Teleflex. Bench test results may not necessarily be indicative of clinical performance. Data on file.

## Turnpike® Catheters

The Turnpike® Catheters are intended to be used to access discrete regions of the coronary and/or peripheral vasculature. They may be used to facilitate placement and exchange of guidewires and to subselectively infuse/deliver diagnostic and therapeutic agents.

The Turnpike® Spiral and Turnpike® Gold catheters are contraindicated for use in vessels with an effective diameter smaller than 1 mm.

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit [teleflex.com](http://teleflex.com).

Teleflex is the home of Arrow®, Deknatel®, Hudson RCI®, LMA®, Pilling®, Rüschi®, and Weck® – trusted brands united by a common sense of purpose.

### Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

### Regional Offices

**United States:** Phone +1 919 544 8000, Toll Free 866 246 6990, [cs@teleflex.com](mailto:cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**Latin America:** Phone +1 919 433 4999, [la.cs@teleflex.com](mailto:la.cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**International:** Phone +353 (0)9 06 46 08 00, [orders.intl@teleflex.com](mailto:orders.intl@teleflex.com), Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

**Australia** 1300 360 226

**Austria** +43 (0)1 402 47 72

**Belgium** +32 (0)2 333 24 60

**Canada** +1 (0)800 387 9699

**China (Shanghai)** +86 (0)21 6163 0965

**China (Beijing)** +86 (0)10 6418 5699

**Czech Republic** +420 (0)495 759 111

**France** +33 (0)5 62 18 79 40

**Germany** +49 (0)7151 406 0

**Greece** +30 210 67 77 717

**India** +91 (0)44 2836 5040

**Italy** +39 0362 58911

**Japan** +81 (0)3 6632 3600

**Korea** +82 2 536 7550

**Mexico** +52 55 5002 3500

**Netherlands** +31 (0)88 00 215 00

**New Zealand** 0800 601 100

**Poland** +48 22 4624032

**Portugal** +351 22 541 90 85

**Singapore** (SEA non-direct sales countries) +65 6439 3000

**Slovak Republic** +421 (0)3377 254 28

**South Africa** +27 (0)11 807 4887

**Spain** +34 918 300 451

**Switzerland** +41 (0)31 818 40 90

**United Kingdom** +44 (0)1494 53 27 61

**For more information, please visit [teleflex.com](http://teleflex.com).**

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions. CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

Teleflex, the Teleflex logo, Turnpike, Arrow, Deknatel, Hudson RCI, LMA, Pilling, Rüschi and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Other names may be trademarks of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 10/2018.

© 2018 Teleflex Incorporated. All rights reserved. MC-004362 Rev 0.1